
GYM[®] *form*

AB SWIFT


- EN** Instruction manual
- ES** Manual de instrucciones
- FR** Mode d'emploi
- DE** Gebrauchsanweisung
- IT** Istruzioni per l'uso
- PT** Manual de instruções
- NL** Bedieningsanleitung

Product description / Descripción del producto / Description du produit / Beschreibung des Produkts
 Descrizione del prodotto / Productomschrijving / Descrição do produto

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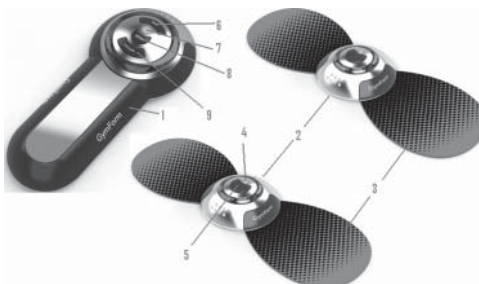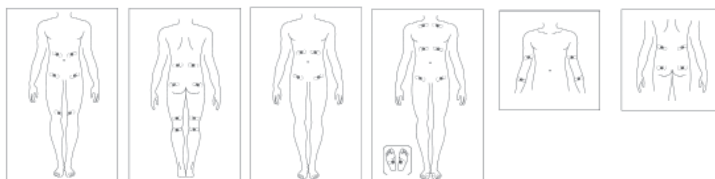Congratulations on your purchase of GYMFORM® AB SWIFT. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your GYMFORM® AB SWIFT. It is essential that anyone who intends to use this device reads and understands the instructions provided. An informed user is a safe user. Thank you for purchasing GYMFORM® AB SWIFT. We hope you enjoy your sessions.

• WARNING

Before embarking on any exercise programme, including the GYMFORM® AB SWIFT, consult your general practitioner and your physiotherapist. Ensure that you read the instructions, warnings and precautions carefully before using your GYMFORM® AB SWIFT.

GYMFORM® AB SWIFT is designed for healthy adults. However, the device should be used with care. Read the following guidelines and the instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.
- Using the GYMFORM® AB SWIFT while connected to high-frequency surgical equipment may cause burns at the contact points of the GYMFORM® AB SWIFT and may result in damage to the product.
- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of the GYMFORM® AB SWIFT.
- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually increase the intensity.
- Do not use on children or animals.
- Do not use while driving or sleeping or if you feel drowsy.
- Do not use in the bathroom or in other places with a high amount of moisture or humidity.
- This device may be used by children aged 8 and over, by people with reduced physical, sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails. Children should not be allowed to play with the device. The cleaning and maintenance to be performed by the user should not be done by children without supervision.
- AVOID using the Gymform® AB SWIFT near devices that emit electromagnetic waves, such as telecommunication antennas, microwave ovens and X-ray equipment. It may interfere with these devices.
- How long the electrodes last will depend on the kind of skin, intensity level used and environmental and storage conditions. For best results, we recommend replacing them after approximately 30 sessions. Spare electrodes are available at the point of sale where you purchased the product. If you notice the electrodes lose power too soon, moisten them with a damp cloth.

DO NOT use GYMFORM® AB SWIFT:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.
- If you are, or think you might be, pregnant.
- If you have cancer, epilepsy or are under medical supervision for a cognitive disorder.

- If you are breastfeeding.
- If you have a high fever.

Women: before using your GYMFORM® AB SWIFT, you should wait:

- At least six months after the birth of your child (after consulting your doctor).
- A month after having an intrauterine device (IUD) fitted.
- At least three months after having a caesarean section (after consulting your doctor).
- Until after the first days of your menstrual period.

Before using GYMFORM® AB SWIFT, you should obtain permission from your doctor or physiotherapist:

- If you suffer from back problems. In this case, the device should only be used at low intensity levels.
- If you suffer from a serious disease or injury not mentioned in this section.
- If you were operated on recently.
- If you are diabetic and take insulin.
- Do not place one receiver on your chest and another on your back in such a way that electricity passes through your heart. Do not place the receiver on genitals, scars, moles, varicose veins or the head.
- Do not place the receiver on recent scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity.
- Do not place the receiver on areas with injuries or limited mobility (e.g. fractures or sprains).
- Do not place the receiver directly on metal implants. It should be placed on the closest muscle.

Possible Adverse Reactions

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne.

If any of these symptoms occur, switch off the device and consult your doctor. If you have never exercised before or have not exercised for a long time: Consult your general practitioner.

Precautions:

- Do not use GYMFORM® AB SWIFT if your body is in contact with water.
- Do not use GYMFORM® AB SWIFT in combination with any other device that transmits electrical current to your body or with any other machine.
- Do not connect the metal snap fasteners or electrodes to any other object.
- Do not modify the product.
- Use only the electrodes supplied by Industex.
- Do not attach two self-adhesive electrodes together.
- If any of the electrodes is not attached properly, you will notice an itching sensation.

SPECIAL FEATURES

1. Compact, wireless easy-to-use design.
2. Working range of up to 2 metres with practically no interference.
3. Thanks to the remote control (radiofrequency technology), the receivers are easy to control regardless of where they are applied on the body.
4. Synchronization of several independent receivers without interference.
5. Synchronization of up to 10 independent receivers without interference. If you acquire additional Gymform AB Swift receivers, you can use a total of 10 at the same time and activate them with the remote control. Make

sure there is a sightline between the remote control and the receivers. If the signal does not reach the receiver, try again by ensuring there is a sightline between the remote control and the receiver.

6. 10 output intensity levels to choose from.
7. 6 different power programmes to choose from.
8. The load-detection function automatically shuts off after 5 seconds without contact with the skin.
9. Made with surface-mount technology to offer better quality and effectiveness.

CONTENTS:

- 2 GYMFORM® AB SWIFT receivers
- 2 Self-adhesive electrodes
- 1 GYMFORM® AB SWIFT remote control
- 1 Instruction manual
- 1 Travel bag

PRODUCT DESCRIPTION:

1. Remote control.
2. Gymform® AB SWIFT receivers.
3. Self-adhesive electrodes.
4. Receiver indicator light.
5. Receiver ON/OFF button. 
6. Button to turn on and increase intensity.
7. Remote control indicator light.
8. Programme button.
9. Button to turn off and reduce intensity.

BEFORE FIRST USE

1. To use the GYMFORM® AB SWIFT, 1 3V CR-2032 battery is required for each receiver (2) and for the remote control (1) (not included). Slide off the cover and insert the battery with the positive end (+) facing up. Close the battery compartment cover.
2. Attach each receiver (2) to each self-adhesive electrode (3) using the 2 snap fasteners and make sure it is properly attached.
3. Remove the protective plastic backing from the self-adhesive electrode (3).
4. Place the self-adhesive electrode (3) on the part of the body to be treated. Figure A provides examples of the areas that can be treated at the same time.
5. Press the receiver ON/OFF button (5) to turn on the receiver. The receiver indicator light (4) will go on to show it is ready to be used.

INSTRUCTIONS

Once you have placed the self-adhesive electrodes (3) on the areas to be treated, you have to turn on the receivers (2) by pressing the button (5) with the symbol , which means ON/OFF. Once you have done this, you can choose from 6 different programmes and 10 intensity levels:

- Button to turn on and increase intensity (6): Press this button to turn on the receivers and increase the stimulation intensity. Each time you press the button, the intensity level will be increased (1-10). There are 10 levels in total.
- Button to turn off and reduce intensity (9): Press this button to turn off the receivers and reduce the stimulation intensity. Each time you press the button, the intensity level will be reduced (1-10). There are 10 levels in total. The receivers can be turned off by pressing the button to the lowest intensity level.
- There are 10 stimulation intensity levels. Level 10 is the one that provides the most stimulation.
- You can change programmes at any time by pressing the

Programme button (8). Each time you press it, a new programme will begin. • There are 6 GYMFORM® AB SWIFT programmes.

Programme 1	Rapid impulses: 18 rapid, consecutive impulses/contractions every 5 seconds/1-second pause between series.
Programme 2	Massage: 16 medium-deep contractions every 10 seconds/1-second pause between series.
Programme 3	Combination of 3 different kinds of contractions every 15 seconds/1-second pause between series. 1) Slow-deep contraction 2) Fast impulse 3) Slow-deep contraction
Programme 4	Combination of 2 different contractions every 2 seconds continuously: 1) Massage: 2 seconds 2) Rapid impulses: 2 seconds
Programme 5	Combination of 4 different contractions every 8 seconds continuously: 1) Rapid impulses 2) Deep contraction 3) Mild contraction 4) Deep contraction
Programme 6	Combination of 8 different contractions every 45 seconds continuously: 1) Slow impulses 2) Deep massage 3) Slow-fast-slow impulses 4) Deep massage 5) Deep impulses 6) Deep massage 7) Mild massage 8) Deep massage

- The maximum recommended time for massaging each part of the body is 10 minutes.
- After using the product, turn off the receivers by pressing the ON/OFF button (5) with the symbol 

Load-detection function

- The receiver indicator light (4) will automatically go out if the device is not in use for 60 seconds. This function prevents the battery from being used if you forget to turn off the receivers.

Automatic timer function

- The device is default programmed to run for 10 minutes. It will automatically shut off when the timer reaches the programmed time.
- When the device is in use (with a stimulation intensity of from 1 to 10), the automatic timer function will start. When the device is on, if the programme is changed while it is in use, the time will maintain the same countdown and the device will automatically shut off when the preset time ends.

TECHNICAL SPECIFICATIONS:

Power Supply: 3 3V CR-2032 batteries

Battery:	3V CR-2032
Maximum voltage:	120V (max.)
Output current:	70 mA (RMS)
Output intensity level:	1-10 (10 levels)

Output programmes:	1-2-3-4-5-6												
Automatic shutoff time:	10 minutes or 1 minute when device is inactive												
Automatic load-free shutoff time:	5 seconds												
Timer:	10 minutes												
Working range:	Distance of 1.5 to 2 m												
Impulse frequency:	<table border="1"> <tr> <td>1</td> <td>5</td> </tr> <tr> <td>2</td> <td>45</td> </tr> <tr> <td>3</td> <td>5-15</td> </tr> <tr> <td>4</td> <td>45</td> </tr> <tr> <td>5</td> <td>4-15</td> </tr> <tr> <td>6</td> <td>4-15</td> </tr> </table>	1	5	2	45	3	5-15	4	45	5	4-15	6	4-15
1	5												
2	45												
3	5-15												
4	45												
5	4-15												
6	4-15												
Impulse width:	200 µs												
Frequency Rx/Tx	2,425 GHz												

Values measured with a load of 10 KΩ.

Temperature range: -45°C to -80°C
 30%-75% RH (without condensation)
 700 -1060 hPa
 Water protection: IPX0.
 Continuously operating unit.
 Pad dimensions: 150 x 50 mm
 Remote control dimensions: 90 mm (length) x 45 mm (width) x 15 mm (height)
 Weight of remote control (without battery): 15 g
 Dimensions of receiver: diameter 45 mm; height: 15 mm
 Weight of receiver (without battery): 10 gr.

This symbol indicates that the device is protected against electrostatic discharge by BF-type double insulation.

ISO 7000-1641 User Guide

CLEANING AND MAINTENANCE:

DO NOT WASH THE Gymform® AB Swift UNIT: WIPE WITH A DRY OR DAMP CLOTH.

Remove the batteries from your GYMFORM® AB SWIFT if you do not plan to use it for an extended period of time.

MAINTENANCE AND STERILIZATION OF ADHESIVE ELECTRODES

- Do not immerse in water.
- Keep in a cool, dry place at room temperature.
- Do not remove the insulating label.
- After using the device, always replace the protective plastic backings on the self-adhesive electrodes.
- Store in the bag supplied.
- Sterilize by wiping with wet fingers.
- Do not remove the insulating label.

PRECAUTIONS WITH BATTERIES

With use of approximately 10 minutes a day, the battery should last for 15 to 20 uses.
 Ensure the batteries are properly installed before turning the unit on.

Turn off the unit to change the batteries.

Remove the batteries if you are not going to use the device for an extended period of time.

Take the necessary precautions when disposing of the batteries. Do not expose them to fire.

CERTIFICATES AND DIRECTIVES

This product complies with Low Voltage Directive 2006/95/EC, Electromagnetic Compatibility Directive 2004/108/EC, RoHS Directive 2002/95/EC, and R&TTE Directive 1999/5/EC.

DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations. Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

The battery should be disposed of separately. To protect the environment, all batteries, including rechargeable ones, should not be disposed of as you would other household waste. Take them to a disposal centre and ensure it complies with applicable legislation on battery disposal.

Quality Warranty:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia only: This product has the benefit of certain Consumer Guarantees prescribed by the Australian Consumer Law which provide protection for consumers. There is no separate warranty as referred to above.

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ISL GmbH
Best Direct International PTY LTD
ISL (Italy) S.r.l.
Best Direct International Ltd,