# beurer

# HK 25, 35 40, 63

HK 25\*

HK 35\*

HK 40

HK 63\*

- (D)** Heizkissen  
Gebrauchsanleitung
- (GB)** Heating pad  
Instruction for Use
- (F)** Coussin chauffant  
Mode d'emploi
- (E)** Almohadas eléctricas  
Instrucciones para el uso

- (I)** Cuscino riscaldante  
Istruzioni per l'uso
- (TR)** Isıtılmalı yastık  
Kullanma Talimatı
- (RUS)** Электрическая грелка  
Инструкция по применению
- (PL)** Poduszka elektryczna  
Instrukcja obsługi

Beurer GmbH  
Söflinger Str. 218 • 89077 Ulm, Germany  
Tel.: +49 (0)731 / 39 89-144 • Fax: +49 (0)731 / 39 89-255  
www.beurer.de • Mail: kd@beurer.de

**Zeichenerklärung Etikett • Explanation of symbols • Légende des symboles**  
**• Explicación de los símbolos • Spiegazione dei simboli Etichetta • Etiket işaret tarifi**  
**• Объяснение условных знаков на этикетке • objaśnienie oznaczeń etykiet**

Anweisungen lesen!  
 Read instructions!  
 Lire les consignes!  
 ¡Lea las instrucciones!

Leggere le istruzioni!  
 Talimatları okuyun!  
 Прочеть инструкции!  
 Przeczytaj instrukcje!

Nicht gefaltet oder zusammengeschoben  
 gebrauchen!  
 Do not use when folded or creased!  
 Ne pas utiliser plié ou comprimé!  
 No utilizar estando plegado o arrugado!  
 Non utilizzare ripiegato né accartocciato!

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**Dear Customer,**

We are happy that you have decided on a product from our range. Our name stands for high-quality and thoroughly tested products from the areas of warmth, blood pressure, body temperature, weight, gentle therapy, massage and air.

Please read these instructions for use carefully and follow the directions.

Best regards,  
The Beurer Team

# 1. Important safety information – read carefully and keep for later use

Failure to observe the information provided below can lead to personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information not only serves to protect your health and the health of others, it also protects the product from damage. Therefore, observe this safety information and pass on these instructions together with the product.

- Only use this heating pad for the purpose described in these instructions.
- This heating pad is not intended for use in hospitals.
- This heating pad is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety.
- This heating pad should not be used with persons who are helpless, small children or persons insensitive to heat (e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol).
- Excessively long use can cause the skin to be burned.
- Caution! Never fall asleep while the heating pad is switched on.
- Do not apply the heating pad to parts of the body which are inflamed, injured or swollen. If in doubt, seek medical advice before use.
- The electric and magnetic fields emitted by this heating pad can interfere with the operation of pacemakers. However, they are far below the permissible limits: electric field strength: 5,000 V/m max., magnetic field strength: 80 A/m max., magnetic flux density: 0.1 millitesla max. Therefore, please consult your doctor and the manufacturer of your pacemaker before using this heating pad.
- Heating pad
  - only connect to the voltage specified on the appliance,
  - HK 25, 35, 63: only use in conjunction with switch types specified on the heating pad,
  - never operate without supervision,
  - never switch on when folded or compressed,
  - never wedge between objects,
  - never make any sharp creases,
  - never use with animals,
  - HK 25, 35, 63: never use when damp,
  - the switch and mains cable should not be exposed to moisture.
- When in use, the electronic components in the switch of the heating pad make the switch slightly warm. The switch should therefore never be covered or placed on the heating pad when it is being operated.

- Do not pull, twist or make any sharp kinks in the cables.
- Do not stick pins or pointed objects into this heating pad.
- This appliance may not be used by children unless the switch has been preset by a parent or a supervisory person, or the child has been sufficiently instructed in how to use this heating pad safely.
- Children should be supervised to ensure that they do not use the heating pad as a toy.
- This heating pad must be checked frequently for signs of wear or damage. If such signs are present, or if the cables are damaged, or if this heating pad has been used incorrectly, it must be taken to either the manufacturer or the dealer before further use.
- Repairs may only be carried out by specially trained personnel or at an authorised repair workshop, as special tools are required. Incorrect repairs may cause considerable danger to the user.
- Should you still have questions on using our appliances, please contact our customer service.

## 2. Proper Use

This heating pad is designed for warming the human body. It allows you to apply targeted warmth. Heat promotes the circulation and relaxes the muscles. This heating pad is not intended for use in hospitals or for commercial applications. In particular, it may not be used on infants, small children, persons insensitive to heat, helpless persons or animals.

To avoid soiling or damaging the heating pad, we recommend using a fabric cover, particularly when applying essential oils and ointments.

When used for the first time, the heating pad may start to smell of plastic but this will soon go away.


## 3. Operation

### 3.1 Safety

This heating pad has a SAFETY SYSTEM. This electronic sensor technology prevents the entire surface of the heating pad from overheating by automatically switching off in the event of a fault. If the safety system switches off the heating pad automatically as a result of a fault, the temperature display on the switch is no longer illuminated even when the heating pad is switched on. Please note that, after a fault has occurred, the back and neck heating pad can no longer be used for safety reasons and must be sent to the service address indicated.

### 3.2 Putting into operation HK 25, 35, 63

First check that the switch is connected to the heating pad at the plug-in connection. Plug the power cord into the mains outlet.
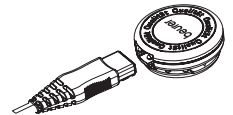
### 3.3 Selecting temperature

The fastest way to heat up the heating pad is to first set the switch to the highest setting. A lower setting can be selected later if required.

After switching on the heating pad, the setting starts flashing briefly, while the heating pad conducts a self-test. Then the setting remains illuminated continuously.

### 3.4 Temperature settings

- Level:
- 0 Off
  - 1 Minimum heat
  - 2 Medium heat
  - 3 Maximum heat

### 3.5 Automatic switch-off function

The heating pad switches off automatically after approx. 90 minutes. The illuminated temperature setting display then starts to flash. To switch on the heating pad again, first set to '0' and then after approximately 5 seconds select the desired temperature setting. Switch off the heating pad (0) or pull the plug out of the mains socket if you are not going to use the heating pad after the automatic switch-off.

### 3.6 Fast heating HK 35, 40, 63

These heating pads have a fast heating function which allows them to heat up quickly within the first 10 minutes.

### 3.7 Additional Instructions for HK 40 Moist heating pad

The HK 40 is also suitable for moist applications. For this purpose it is equipped with an additional moisture-storing insert of fleece material. The heat transfer to the body is considerably improved during moist use. Compared to a conventional, dry use of the heating pad, the moist heat is generally perceived as more pleasant.

**For safety reasons, the textile covering and the insert may only be moisturised while the heating pad is disconnected from the power supply (with the power plug removed from the electrical outlet).**

To use when damp, sprinkle or spray the textile covering and the insert with water.

**Before switching on the heating pad, always make sure that the textile covering and the insert no longer drip, even when subjected to a higher press-on pressure! The switch must never come into contact with water. Never allow water from an excessively damp covering to run along the cable into the switch.**

The plastic casing of the heating pad is watertight. Please make sure that the sleeve is never damaged and do not use the heating pad if signs of damage or wear appear.

When using the (damp) heating pad, no essential oils or creams may be used.

These can result in damage to the plastic casing under certain circumstances.

You can also use the heating pad dry (without moistening it).

### 3.8 Additional information for HK 63 (CE 0344)

Medical product: This heating pad has 6 integrated permanent magnets, which generate a homogeneous magnetic field. Under certain conditions, static magnetic fields can have a favorable effect on the organism and relieve pain. This heating pad can provide relief particularly in rheumatic complaints. However, the heat pad must not be used in inflammatory diseases.

## 4. Cleaning and Care

### 4.1 HK 25, 35, 63

- **Before cleaning the heating pad, always unplug it from the mains.**
- The textile cover can be cleaned according to the symbols on the label, and for this purpose must first be removed.
- Please note that the actual heating pad may not be machine-washed, spun or machine dried, chemically cleaned, wrung out, pressed or ironed.
- Smaller stains can be removed using a cloth or moist sponge together with a little liquid detergent for delicate fabrics if necessary.
- **Only use the heating pad again once the heating pad and cover are completely dry.**

### 4.2 HK 40

Moist heat can promote the growth of germs. Therefore, make sure after each damp use that the heating pad dries again completely and the covering and insert are washed at frequent intervals. The insert can be washed in the washing machine in hot water (approx. 140 °F). The textile covering can be cleaned in accordance with the symbols on the label. Remove it for this purpose. Also clean the plastic casing of the heating pad from time to time with a damp cloth.

## 5. Storage

If you are not going to use the heating pad for some time, we recommend storing it in its original packaging in a dry place without weighing it down.

Allow the heating pad to cool down beforehand.

## 6. Disposal

Please dispose of the blanket in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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