FOR A GOOD REASON

GRUNDIG

**MULTI-FUNCTION
BODY ANALYSER
SCALES**

PS 5110

de

en

fr

it

pt

es

sv


nl

pl

cs

tr

el

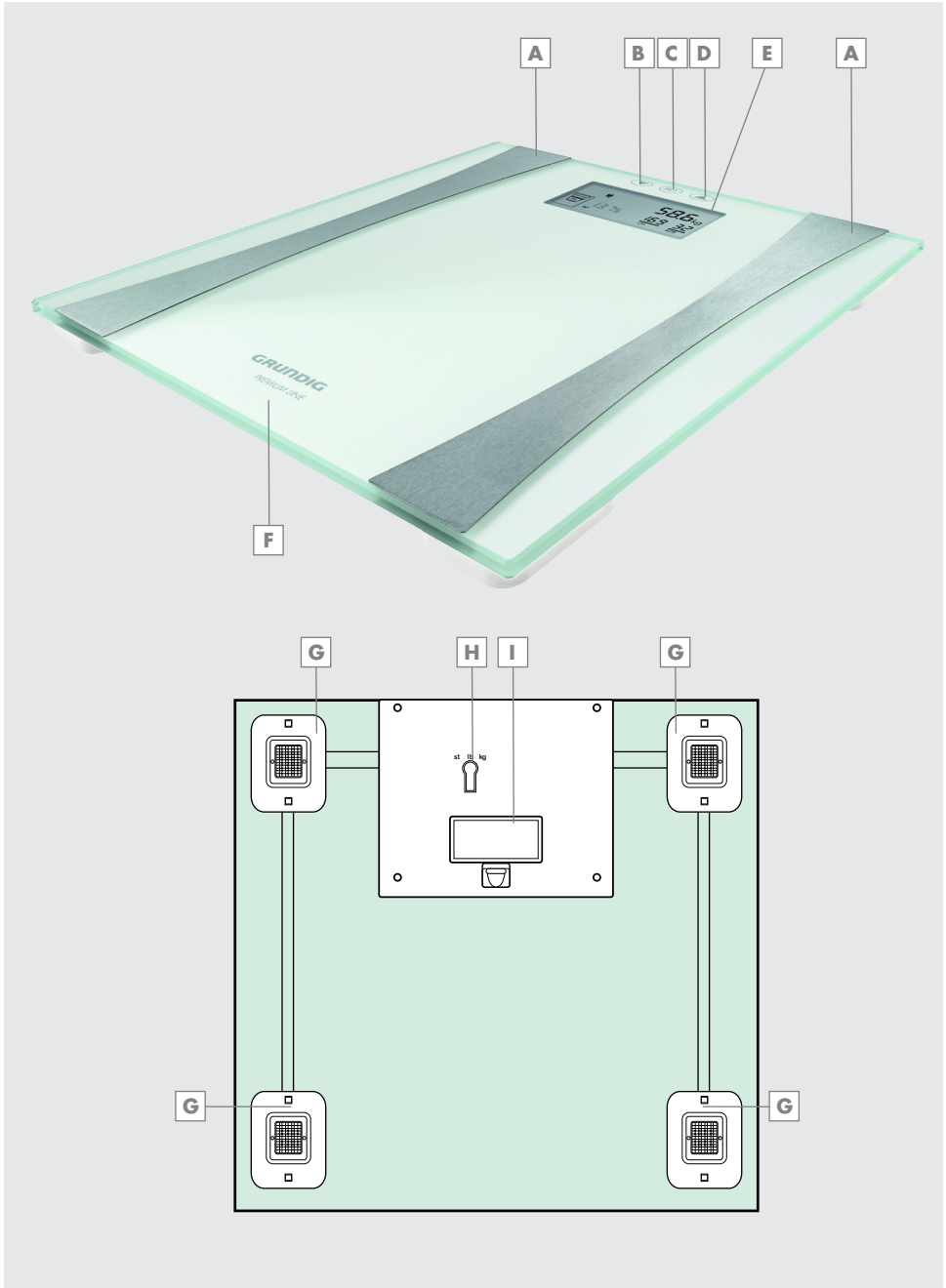DEUTSCH 04-12

ENGLISH 13-20

FRANÇAIS 21-28

ITALIANO 29-36

PORTUGUÊS 37-44

ESPAÑOL 45-52

SVENSKA 53-60

NEDERLANDS 61-68

POLSKI 69-76

ČESKY 77-84

TÜRKÇE 85-92

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SICHERHEIT

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SAFETY AND SET-UP

When using the scales for the first time, please pay attention to the following:

- This scales is only intended for household use.
- Do not use the scales on a carpet or carpeted flooring.
- Always use the scales on the same smooth and hard floor. This ensures consistent measured values.
- When stepping on the scales, keep a distance of 2 cm to the edge. Please be careful when doing this.
- First centralise your body weight in the middle of the scales before taking one foot off the scale. Here, always make sure that the scales cannot tip over (risk of injury).
- Never step on the scales with wet feet. Danger of slipping!
- Never use the scales outdoors.
- Stop using the scales if defects occur.
- Stop using the scales if the glass surface breaks or has a crack in it.
- Do not place any objects on the glass surface of the scale.
- Never jump on or off the scales.
- If the scales is used in the bathroom, ensure that it does not come into contact with water.
- If you have problems with this device, try resetting it by removing the batteries briefly or replacing them.
- Never immerse the scales into water.

- Do not expose the scales to direct sunlight, high temperatures, or constant high humidity.
- Keep the scales away from children.
- Do not operate the scales if it shows signs of visible damage.

- Do not open the scales under any circumstance. For damage occurring as a result of improper handling, any claims under guarantee are no longer valid.

Caution

- The scales may not be used by the following groups:
Pregnant women and youths under 18 years old should not use the scales for body fat analysis (in such cases, an exact body fat value cannot be determined).
- The scales measures your proportion of body fat using BIA (Bioelectrical Impedance Analysis). To do this, a harmless signal is sent through your body. Thus this scales should not be used by persons with a pacemaker or with any other kind of operatively inserted medical device. When in doubt, consult your physician.
- Your product is designed with complex algorithms. Due to these algorithms, if an external electromagnetic disturbance happens near your device, it may affect your readings. In such cases, first remove the electromagnetic disturbance reason and repeat your measurement. Under extreme electromagnetic disturbance, your device will protect itself and enter into protection mode. To turn back to normal running mode, remove and replace the battery of the product after removing the disturbance.

OVERVIEW

Dear customer,

Congratulations on the purchase of your PS 5110 Multi-function Body Analyser Scales.

Please read the following user guide for this scales thoroughly so that you may enjoy many years of use with this high-quality Grundig product!

A responsible approach!

GRUNDIG focusses on contractually agreed social working conditions with fair wages for both internal employees and suppliers. We also attach great importance to the efficient use of raw materials

with continuous waste reduction of several tonnes of plastic every year. Furthermore, all our accessories are available for at least 5 years. For a future worth living.

For a good reason. Grundig.

Features

In addition to displaying your body weight, this scales also enables a comprehensive body composition analysis.

Depending on the selection of various preset user parameters and the input of individual values, the following values and information are calculated and displayed:

- BMI (Body Mass Index)
- Percentage of body fat - TBF
- Weight interpretation – slim, normal, slightly overweight, overweight
- Percentage of water in body tissue - TBW
- Your ideal weight, based on the data you entered.
- Percentage of muscle weight - TBM
- Percentage of bone weight in comparison to body weight in kg (bone mass)
- A daily calorie intake suggestion

Caution

- The PS 5110 is a product for your well-being. It is not a medical product. No medical recommendations can be made on the basis of its use or of the results displayed.

Brief description

Refer to the diagrams on page 2.

- A** Measuring electrodes (left and right)
- B** Button  (on narrow top side)
- C** Button  (on narrow top side)
- D** Button  (on narrow top side)
- E** LCD screen
- F** Glass surface
- G** Sensors with anti-slip stands
- H** Unit selector switch (st, lb, kg)
- I** Battery compartment

STARTING OPERATION

Power supply

This scales operates on 2 x 1.5 V LR3/AAA batteries (total of 3 V).

- 1 Open the cover of the battery compartment .
- 2 Insert the 1.5 V LR3/AAA batteries. When doing so, pay attention to the polarity markings in the bottom of the battery compartment.
- 3 Close the cover of the battery compartment .

Notes

- Do not expose the batteries to extreme heat e.g. direct sunlight, heaters, fire, etc.
- When the batteries are weak, » L o « will appear in the display. Replace the batteries immediately.
- Remove the batteries when they are flat or when you know that the device will not be used for a longer period of time. No liability will be assumed for damage incurred as a result of leaking batteries.

Environmental note

- Batteries, including those which do not contain heavy metal, should not be disposed of with household waste. Please dispose of used batteries in an environmentally sound manner. Find out about the legal regulations which apply in your area.

Operating as a scales

Notes

- You do not need to enter any data for measuring your body weight.
 - The scales measures your weight in increments of 100 grams.
 - It uses a very precise 4 point measurement.
- 1 On the bottom of the scales, use the switch  to select between Imperial (st/lb/in) and European measurements (kg/cm).
 - 2 Position the scales on a smooth, flat surface.
 - 3 Activate by tapping your foot once on the scales in the bottom right corner of the glass surface .
 - Display: » 888 « appears on the LCD screen .

- 4 As soon as the display shows » 0.0 kg « carefully step on the scales.

- 5 Stand still on the glass surface until the weight displayed flashes twice and then no longer changes.
 - Your body weight is shown in the display.
 - After about 5 seconds, the scales automatically switches off.

Note

- If the » 0-L d « display appears during weighing, the scales is overloaded. After about 5 seconds, the scales automatically switches off.

OPERATION

Operating as a body analyser

The measurement and data is transmitted through your bare feet, which are positioned on the measuring electrodes **A**. Therefore, ensure that both feet are placed on the scales properly.

Entering your data

When you use the body composition analysis function for the first time, select the following data and then enter the following information: your user parameter (from 01 to 12), your gender, your height (75 - 215 cm), your age (10-99) and activity level (select from A 1 to A 4).

Notes

- After entering the initial data, you will not have to press the **C** (SET) button again the next time you use the scales. Simply select the corresponding parameter (01 to 12) using the **B** (▲) and **D** (▼) buttons.
- During data entry, if no button is pressed for more than 20 seconds, the scales automatically switches off. You must then start from the beginning.

Entering parameters

- 1 Position the scales on a smooth, flat surface and activate it. To do this, tap your foot once on the scale in the bottom right corner of the glass surface **F**.
 - Display: »**BBB**« appears on the LCD screen **E**.

- 2 Press the **C** (SET) button as soon as »**00** kg« appears in the display.
 - Display: »**0** t« flashes.

Note

The selection 01 is described below as an example.

- 3 Select the desired user parameter (01 to 12) using the **B** (▼) and **D** (▲) buttons and confirm your selection with the **C** (SET) button.
 - The symbol  flashes.

- 4 Confirm the selection of a male user with the **C** (SET) button.
 - Display: The height flashes.

- 5 Set the height with the **B** (▼) and **D** (▲) buttons and confirm your selection with the **C** (SET) button.
 - Display: The age flashes.

OPERATION

- 6 Set your age with the **B** and **D** buttons and confirm your selection with the **C** button.
 - Display: Activity level A1 flashes.

- 7 Set your personal activity level with the **B** and **D** buttons:

A 1	Unfit person used to little movement
A 2	Person who mostly works sitting down
A 3	Person used to slight physical exertion
A 4	Person used to hard physical work, such as dancers, athletes, etc.

- 8 Confirm the selection with the **C** button.

Performing the analysis

After entering your parameters, you can use the body composition analysis function of the scales. This process is described below using 01 as the user parameter:

- 1 Position the scales on a smooth, flat surface, and activate it. To do this, lightly tap your foot once on the scales in the bottom right corner of the glass surface .
 - Display: »888« appears on the LCD screen .

- 2 As soon as »0.0 kg« appears in the display, select the desired user parameter (01 to 12) using the **B** and **D** buttons.


Notes

- The values saved for this person are displayed one after another, and then the scales switches to weighing mode.
 - Display: »0.0 kg« appears in the display .
- If no data has been saved, the scales immediately switches to weighing mode.
 - Display: » 0.0 kg « and » « appear in the display .

- 3 With bare feet carefully step onto the glass surface of the scales. Make sure that your feet are positioned on both measuring electrodes .

- The scales displays your weight and the » « symbol rotates in a circle twice.

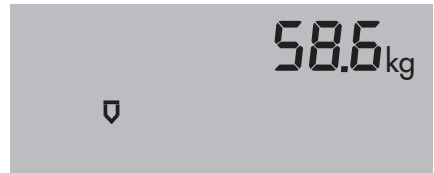OPERATION


Notes

The following values and information are shown one after another in the display **E**:


- BMI (Body Mass Index) :


- The percentage of body fat - TBF, as well as a weight interpretation - slim, normal, slightly overweight, overweight.

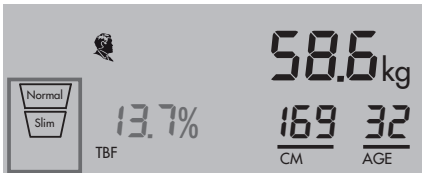- If the body fat percentage is unrealistic, »Err« will be shown in the display and the scales will automatically switch off.
- The percentage of water in your body tissue - TBW, as well as your ideal weight, based on the data you entered.

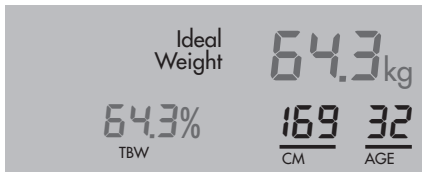- The percentage of muscle weight - TBM, the percentage of bone weight in comparison to body weight in kg (bone mass), as well as a suggestion for daily calorie intake (Kcal).

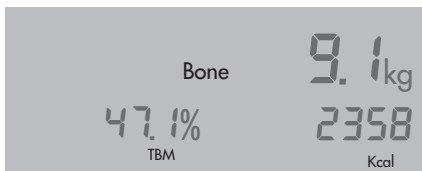- The values and information provided are displayed twice and remain visible in the display for approx. 6 seconds each time. After that, the scales automatically switches off.

Notes

- Separate values are saved for each user (01-12).
- The data is deleted when you remove the batteries.

NOTES ON MEASUREMENTS

Total body water

Water is an essential part of the body. The amount of water in the body may vary between 35% and 70%. It helps all the essential biochemical reactions necessary for normal bodily functions. Water transports waste products out of the cells and discharges them in sweat and urine.

Water also maintains a constant body temperature. The amount of water in the body depends on the proportion of water loss, for example due to heat or physical activity and sweating. Therefore, it is very important to check the proportion of water in the body and ensure that it remains constant.

Notes

- The following tables give you a reference which can help you to judge your body fat and body water yourself. There is no standard reference data for muscle mass or bone weight.
- It is important to see and follow the development of your body fat, body water, muscle mass and bone weight over a period rather than concentrate on the body weight only.

- Make sure that your bare feet are well positioned on the stainless steel electrodes. Extremely dry skin or thick callus could influence the analysis results. This may prevent the scales from functioning properly and the error message »Err« will be displayed.

Body fat composition

Checking the body fat composition is particularly important. A high fat mass can lead to high blood pressure, high cholesterol, heart disease and other illnesses. The American College of Sports Medicine have made the following recommendations for body fat mass:

Notes

- Measure body weight, body fat percentage and water proportion at the same time every day, preferably in the evening.
- Use the scales when you are undressed if possible. This will give you more exact results.
- The results are not suitable for diagnosing medical problems. If in doubt, seek medical advice.
- Your feet should always be clean and dry.

Age	Body fat (women)				Body fat (men)			
	Slim	Normal	Slightly over-weight	Over-weight	Slim	Normal	Slightly over-weight	Over-weight
19 - 39	<20%	20-32%	33-38%	>38%	<12%	12-20%	21-25%	>25%
40 - 55	<23%	23-35%	36-41%	>41%	<13%	13-21%	22-26%	>26%
56 - 85	<24%	24-36%	37-42%	>42%	<14%	14-22%	23-27%	>27%

Age	Body water (women)			Body water (men)		
	Low	Normal	High	Low	Normal	High
18 - 39	47	47-57	57	53	53-67	67
40 - 55	42	42-52	52	47	47-61	61
56 - 85	37	37-47	47	42	42-56	56

Notes on operating and cleaning

- In order to preserve the batteries, the scales is equipped with an Auto-off function. It automatically switches off when operation has been interrupted for more than 20 seconds.
- If you would like to achieve accurate measuring results, do not use the scales directly after consuming coffee, beer, wine or other alcoholic beverages. The same applies for extreme physical exertion. These factors reduce the percentage of water in body tissue (dehydration). In such cases, wait a few hours before using the scales.
- Clean the scales with a damp cloth. Do not use any cleaning agents. Do not submerge the device in water.

Technical data

This product complies with the European Directives 2004/108/EC.

Power supply

2 x 1.5V LR3/AAA batteries, 3 V

Maximum weight load

150 kg

Minimum weight load

2.0 kg

Ambient temperature

20 °C ± 15 °C

Technical and design modifications reserved.

Environmental note

This product has been made from high-quality parts and materials which can be re-used and recycled.

Therefore, do not throw the product away with normal household waste at the end of its service life. Take it to a collection point for recycling electrical and electronic devices. This is indicated by this symbol on the product, in the operating manual and on the packaging.

Please find out about collection points operated by your local authority.

Help protect the environment by recycling used products.

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