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# English

## 1. Introduction

**Thank you for your decision to purchase the Chi Swing Massager.**

New technology and the search for different types of massaging techniques are constantly being researched and developed. We feel this massage effectively counteracts the strenuous rigors of everyday life while helping one to relax.

Enjoy the innovative massage of the Chi Swing massager for a relaxing massage on the neck, shoulders back or thighs. The Chi Swing massager will spoil you with its deep penetrating rotating vibration massage. It will loosen tight muscles and help ease tenseness. Circulation will be increased

Choose between the automatic program or the Manuel mode and five massage intensities. The settings can easily be selected or changed from the hand controller. The various settings are indicated on the hand held controller. We hope that the Chi-Swing Massager will provide you with much pleasure.

We hope you have much pleasure with your Chi-Swing Massager.

## 2. Main Features of the Device

### Chi-Swing Massager

- Rotating, deep penetrating vibration massage
- 5 massage intensity settings
- Selectable between Manual or preprogrammed massage settings
- Can be placed 3 different positions on the back for different massage effect.
- Can also be used on the stomach or upper thigh
- Secure solid hold in each position
- Automatic shut-off after 10 minutes

### Safety precautions

Read the operating instructions carefully!

Keep the operating instructions safe and allow other users access to them as well.

Should you loan the device please include the instructions as well.

To ensure maximum safe operation of the Chi-Swing Massager:

- Before connecting the power supply to the wall socket, confirm that the electric voltage indicated on the specification plate matches your power supply network.
- Pull power plug if the device is not used or during cleaning.
- Do not operate the device if cables or plugs are damaged, if it does not work correctly, or if damaged.
- Keep device and cable away from hot surfaces, heat, moisture and liquids. Do not touch the power plug with wet or humid hands or if you are standing in water.
- Do not reach for a device having fallen into water.
- Do not hold, carry or lift the device at the power supply cable and do not use cables as a handle or for hanging up the device.
- Prevent the device from coming into contact with sharp objects.
- Do not leave the device unattended while in operation.
- Do not let children play with packaging as there is the risk of suffocation.

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- Only use this device for the intended purpose as described in these operating instructions.
- Longer massage periods cause an increased motor heat and thereby a diminished lifespan of the product.

**Use the Chi-Swing Massager indoors only.**

**Do not use the Chi-Swing Massager in damp rooms (e.g. bathroom or prior to taking a shower).**

### Intended application

The Chi-Swing Massager is specifically intended for personal massage of the neck, shoulders, back and upper thigh regions. It is not intended to treat medical conditions and should not replace medical treatment.

**This appliance is only intended for the use in the home and not for commercial or medical purposes. Should you have any medical reservations, please contact your doctor before using the Chi-Swing Massager.**

Do not use the Chi-Swing Massager:


- in case of pacemakers or implants
- in case of a slipped disc or pathological changes of the vertebrae
- in case of bad circulation or headache
- in case of tuberculosis, tumors, dermatitis and phlebitis, thrombosis, bleedings, open and fresh wounds, chapped skin, varicose veins, neuropathies (e. g. sciatica)
- after an operation on the upper part of body, gluteal region, thighs or if you suffer a serious illness
- during pregnancy
- Use is not recommended for children younger than twelve (12).

### 3. Package Contents


Chi-Swing Massager  
12V Power supply  
2 weights  
extension 60 x 5cm  
Instructions D / F / GB / NL / I / E

## 4. Description

### Chi-Swing Massager

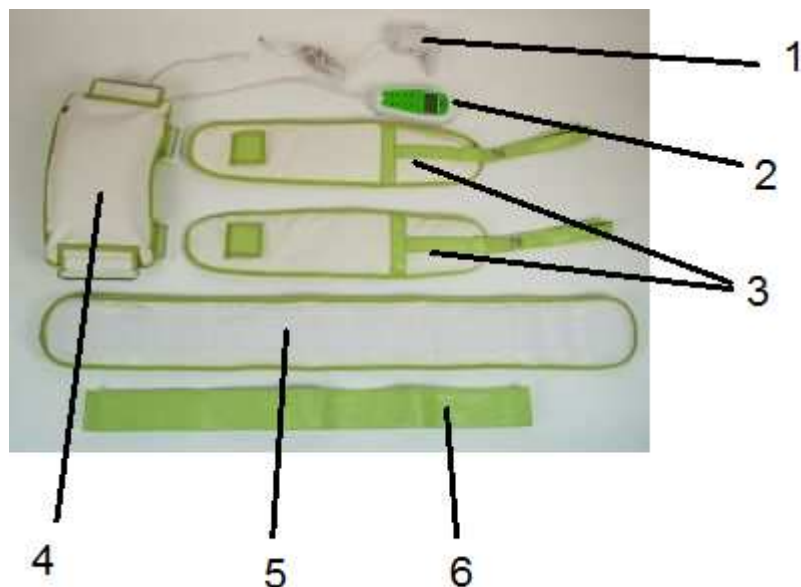1. Mains adaptor
2. Controller
3. Weight belt for Neck/Back massage
4. Chi-Swing Massager
5. Horizontal belt
6. Horizontal belt extension
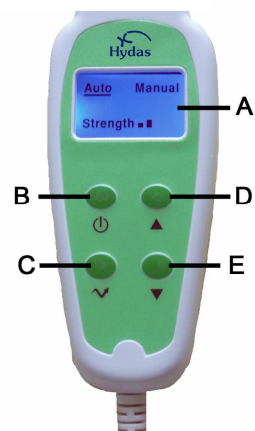
#### Controller

- A – Display  
 B – On / Off button  
 C – Mode selector button  
 D – massage intensity increase  
 E – massage intensity decrease

**Before connecting the massagers plug to the power supply system, make sure the Chi-Swing Massager is switched off. Please use only the adapter included.**

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## Before operating:

- Place the Chi-Swing Massager for the neck, shoulders or back in desired position.

## Alternative:

### Thigh, Stomach, and Lumbar massage

- Place the horizontal belt through the metal rings of the Chi-swing Massager and close
- Place the Chi-swing Massager in desired location, stomach, legs, or hip area

If necessary, the included belt extension gives one the ability to customize their fit for a hip massage.

- Plug circular plug of main adaptor into the provided round slide on receptacle of the Chi-Swing Massager
- Connect the plug of the power supply to the electrical outlet.

## Time of use / breaks


The recommended duration of use is approximately 10 minutes. A break of approximately 1 hour is recommended before further use.

It is not recommended to carry out more than 3 massages/day.

## 5. Application

Turn the Chi-Swing Massager on by pressing the on / off button (B).

The display will light in blue and display the current massage setting as shown in the illustration.

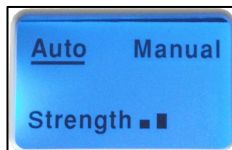**In Program Mode the Chi-Swing massager starts at the second intensity setting of the program.**


Depressing the on / off button (B) again turns the unit off.

The massage mode can be changed manually by pressing the Mode selector button (C) .

By pressing the mode button (C) again returns the unit to the program mode.

### Program Mode

The sequence of the programmed massage mode corresponds to the diagram:

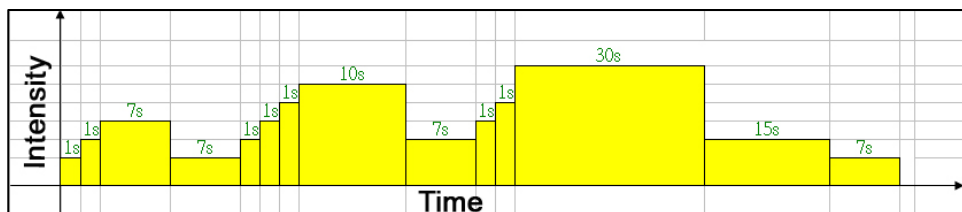The massage intensities vary automatically in Program mode.

The message intensity can also be changed by pressing either button ▲ (D) or ▼ (E).

To increase intensity press button ▲ (D).

To decrease intensity press button ▼ (E).

### Manual Mode

When starting manual mode or switching to manual the massager will begin at the lowest intensity.

To increase message intensity press ▲ (D).

Press ▼ (E) to decrease message intensity.

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During use the mode can be changed at any time by pressing button (C).

In Program Mode the Chi-Swing massager starts at the second intensity setting of the program.

In manual mode, the Chi-Swing massager starts at the first intensity setting.

**Example:** you have selected the fourth massage intensity in the program mode and wish to change to manual mode. Press mode button (C). The Chi-Swing Massager will change to manual mode and begin at the first intensity, by pressing (C) again the massage intensity will increase to the second level.

## Automatic shut-off

The Chi-Swing Massager switches automatically off after 10 minutes in operation. As the recommended massage time is approx. ten minutes, the massager should not be switched on again immediately after shutdown. A rest of about one hour is recommended and three massages per day at the most should not be exceeded.

## 6. Cleaning and Storage

### Cleaning:

Pull power plug after switching off of the Chi-Swing massager and wait until the device has cooled. Only clean Chi-Swing massager and the controller by means of a soft, damp cloth. For this purpose you can apply a mild cleaning agent.

### Storage:

Keep the Chi-Swing massager in a dry, cool place.

It is recommended to store the Chi-Swing massager in its package in order to avoid a contact with sharp edges or objects that could damage the cover.

## 7. Technical Specifications:

|                    |                                  |
|--------------------|----------------------------------|
| Model/Type:        | LH-113R                          |
| Voltage:           | 100-240V ~ 50/60Hz               |
| Mains adaptor:     | 12V                              |
| Plug:              | Euro plug                        |
| Performance:       | 15W                              |
| Material:          | 100% PU (Polyurethan)            |
| Weight:            | 2x á 270g                        |
| Product dimension: | approx. 78 x 15 x 8cm            |
| Weight:            | approx. 2 kg (without packaging) |

## 8. Warranty:

The guarantee lasts for 24 months from the date of purchase for manufacturers' defects. This guarantee is only valid if the instructions in the "instructions for use" have been followed. The guarantee is no longer valid if the device is damaged due

to, misuse, improper use or service by unqualified persons.

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Any questions?

Our hotline is: +49 (69)95406124

For the sake of the environment, we request that you take packaging material to be recycled (collection points).

### Disposal of appliance

Electrical and electronic devices may not be disposed of with domestic waste. Consumers are obliged by law to return electrical and electronic devices at the end of their service lives to the public collecting points set up for this purpose or point of sale. Details to this are defined by the national law of the respective country. This symbol on the product, the instruction manual or the package indicates that a product is subject to these regulations. By recycling, reusing the materials or other forms of utilising old devices, you are making an important contribution to protecting our environment.

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© by Hydas GmbH & Co. KG, Am Hohlen Weg 37, D-34369 Hofgeismar, Germany