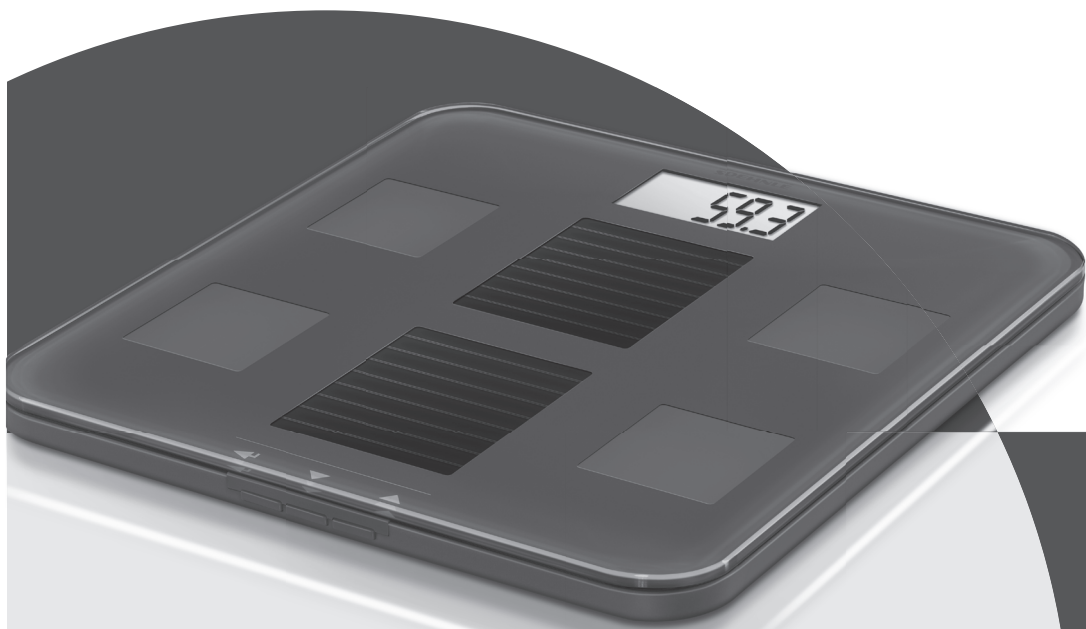
SOEHNLE

SOLAR FIT

BEDIENUNGSANLEITUNG	2
OPERATING INSTRUCTIONS	3
MODE D'EMPLOI	4
ISTRUZIONI PER L'USO	5
GEbruIKSAANWIJZING	6
INSTRUCCIONES DE MANEJO	7

Alter/Age	  					  				
	-	✓	+	++	✓	-	✓	+	++	✓
10-12	< 12	12-23	23-30	> 30	> 60	< 8	8-18	18-24	> 24	> 64
12-18	< 15	15-25	25-33	> 33	> 58,5	< 8	8-18	18-24	> 24	> 63,5
18-30	< 20	20-29	29-36	> 36	> 56	< 8	8-18	18-24	> 24	> 62,5
30-40	< 22	22-31	31-38	> 38	> 53	< 11	11-20	20-26	> 26	> 61
40-50	< 24	24-33	33-40	> 40	> 52	< 13	13-22	22-28	> 28	> 60
50-60	< 26	26-35	35-42	> 42	> 51	< 15	15-24	24-30	> 30	> 59
60 +	< 28	28-37	37-47	> 47	> 50	< 17	17-26	26-34	> 34	> 58


Körperfettanteil, body fat content, taux de graisse, percentuale di grasso, lichaamsvetaandeel, procentaje de la grasa corporal, procentagem de gordura corporal, fedtprocent, andelen kroppsfett, ehon rasvapitoisuutesi, test zsiřtartalmát, zawartość tłuszczu w ciele, podíl tělesného tuku, Содержание жира в организме, vücuttaki yağ oranını, το ποσοστό λίπους, delež tělesne maščobe, udjel tjelesne masnoće.





Körperwasseranteil, body water content, taux d'eau dans le corps, percentuale acqua corporei, lichaamswateraandeel, procentaje del agua corporal, procentagem de água corporal, vandprocent, andelen kroppsvatten, kehon nestepitoisuutesi, test víztartalmát, zawartość wody w ciele, podíl vody v těle, Содержание воды, vücuttaki su oranını, το ποσοστό νερού, delež tělesne vode, udjel tjelesne tekućine.



A


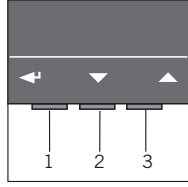
B

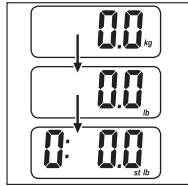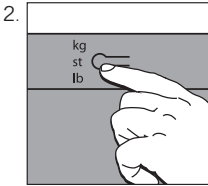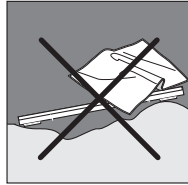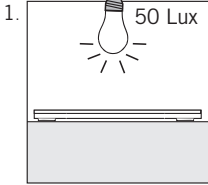
C

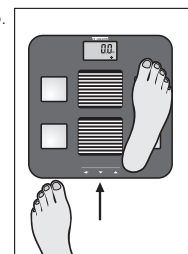
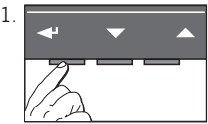
D

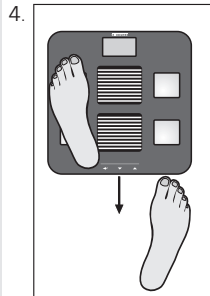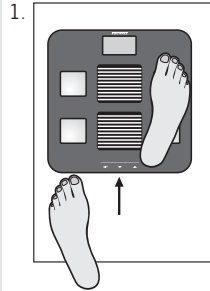E

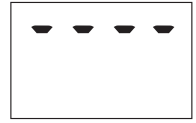F

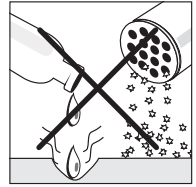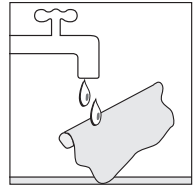A. Bedienelemente

1. Bestätigen (◀)
2. Minus (▼)
3. Plus (▲)

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1. Confirm (◀)
2. Minus (▼)
3. Plus (▲)

B. Preparation

1. For all measurements place scale level on a hard surface (not carpet). Ensure sufficient illumination. Do not cover solar cells.
2. A country-specific option to switch from kg/cm to lb/in or st/in is possible by pushing the key on the bottom of the scale to switch units.

Attention! Slipping hazard in case of wet surface.

C. Data input

Note: Place the body analysis scale on a level surface. Then begin with data entry. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.

The same might occur if the scale is held in the hand prior to the weighting process.

1. Start data entry (push ◀ - key for 1 second)
 2. Select memory slot (P1 ... P8) (▼ = minus, ▲ = plus) and confirm (◀).
 3. Enter body height and age (▼ = minus, ▲ = plus) and confirm (◀).
 4. Select gender using ▼/▲ (♂ = male, ♀ = female) and confirm (◀).
 5. The number of the memory slot is being displayed. Push ◀ key.
- 0.0 appears on the screen.
6. Then immediately step on the scale with **bare feet** for the first body analysis. If you step on the scale with socks or shoes, then only the body weight will be displayed.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data entry must be repeated.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

D. Body analysis

1. Step on the scale with bare feet; stand still for approx. 5 sec.
2. Automatic recognition of persons. Select personal memory slot (▼/▲) if person cannot be automatically recognized.
3. Display weight, body fat content, body water content, and personal memory slot.

Note: Body fat and water are calculated based on the determined weight and the personal data.
4. Step off the scale.

E. Messages

1. Overload: Over 150 kg.

F. Cleaning and Care

Clean only using slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.

G. Technical data

Load capacity x graduation = max. 150 kg x 100 g

Graduation of body fat content: 0.1%

Graduation of body water content: 0.1%

Body height: 100-250 cm

Age: 10-99 years

8 individual memory slots

H. Guarantee

For 3 years as of the date of purchase, Soehnle guarantees the removal of defects due to material or processing errors, free of charge, by way of repair or replacement. In case of a guarantee claim, please return your scale to your seller, together with purchase receipt and guarantee card.

I. EC Conformity

This device complies with the applicable EC Directive 2004/108/EC.

K. Disposal of electric and electronic devices EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

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LEIFHEIT AG