Bedienungsanleitung  
Mode d'emploi  
Istruzioni per l'uso  
Instructions for use  
Instrucciones de uso

**DE – Körperfettwaage**  
**FR – Balance impédancemètre**  
**IT – Bilancia per la misurazione del grasso**  
**EN – Body fat scales**  
**ES – Báscula para medición de grasa corporal**

- 100 g genau, bis 180 kg  
Précision 100 g, jusqu'à 180 kg  
Divisione precisa di 100 g, fino a 180 kg  
100 g accuracy, up to 180 kg  
100 g de precisión, hasta 180 kg
- Präzises 4-Punkte Messsystem  
Système de mesure précis à 4 points  
Precisa tecnologia a 4 sensori  
Precision 4 point measurement system  
Preciso sistema de medición de 4 puntos
- Einfaches Übertragen der Messergebnisse auf Smartphone  
Transmission simple des valeurs mesurées vers le Smartphone  
Facile trasmissione dei dati di misurazione sullo smartphone  
Easy transmission of measurement results to smartphone  
Transmisión sencilla de los resultados de medición al smartphone

**Art. 1862.42**

- 4** Sicherheitshinweise  
Consignes de sécurité  
Direttive di sicurezza  
Safety instructions  
Indicaciones de seguridad

**15**

- Gebrauchen  
Utiliser  
Uso  
Use  
Uso

- 32** Garantie-Hinweis  
Informations de garantie  
Dichiriazione di garanzia  
Warranty information  
Garantía - Nota

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DE

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## Safety instructions

EN

- Children aged 8 or over and people with limited physical, sensory or mental abilities or with and insufficient knowledge and experience should only be allowed to use, clean or maintain the appliance if supervised by someone responsible for their safety or if they have already been shown how to use the appliance safely and have understood the risks involved.
- Children should not be allowed to play with the appliance. Keep the appliance out of reach of children aged under 8.
- Keep packaging (e.g. plastic bags) away from children.
- Never put the appliance on hot surfaces or near open flames. Do not store/expose the appliance to intense heat (from radiators, prolonged sunshine). Only extinguish burning appliances using a fire blanket.
- Do not use attachments not recommended/sold by the manufacturer.
- The appliance is designed for household use and not for industrial operation. Do not operate the appliance outdoors.
- Use appliance in an upright position on a stable surface.
- Regularly check if the appliance is defective. Never switch on damaged appliances – have repairs carried out/obtain replacement from the manufacturer, the manufacturer's service point or from qualified experts. Never open the appliance yourself – risk of injury!
- Do not drop device/knock device against hard objects.
- Never immerse appliance in water/other fluids.
- Do not use in the vicinity of devices with electromagnetic radiation (microwave, telephone, etc.).
- Have the functioning of the device checked every 2 years by a Trisa Electronics service point.
- Insert batteries correctly. Do not mix battery types/new batteries with used batteries. Long-life alkaline batteries are recommended.
- Remove batteries when not in use for a protracted period.

## Safety instructions

EN

-  • Do not use this appliance near the water (bath-tubs, washbasins, etc.) nor expose it to rain or other humidity. Ensure hands are dry when using appliance!
  - If the device falls into water, have it checked out by an electrician before using it again.
-  • Only trained experts should repair electrical appliances. Unskilled repairs can cause considerable danger to the user.
-  • No warranty for any damage can be accepted if the appliance is used for improper purpose, wrongly operated or unprofessionally repaired. In such a case, any warranty claim is void.
  - Technical measuring tolerances are possible as these scales are not calibrated for professional, medical use.
  - Self measurement is for checking only, not for diagnosis or treatment.
  - Do not use the scales when ill or while taking medication that can affect body water percentage.
  - Scales not suitable with: heart pace maker, electromedical equipment, edeme, fever, osteoporosis, dialysis, heart complaints, pregnancy, athletes who train for more than 10 hrs. per week.
  - Do not stand on the scales with wet feet – danger of slipping!

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No household waste. Must be disposed of in accordance with local regulations.

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In accordance with the European guidelines for safety and EMC.

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We reserve the right to make changes in design, equipment, technical data as well as errors.

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Please read all the information provided in these instructions. Keep the instructions in a safe place and pass them on to further users. The appliance must only be used for the intended purpose according to these instructions. Observe the notes on safety regulations.

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Based on age, height and sex, these scales determine within seconds, body weight, body water and body muscle using bio-electric impedance analysis. The measurement current is harmless. The scales are for private use only. Seek advice from your doctor if you have health problems.

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For optimum comparability, always weigh yourself at the same time of day and on the same surface.

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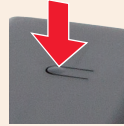
Switch weight unit if required (remains stored):

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Scales only display weight. Body fat, body water, muscle mass and bone mass are displayed in the app.

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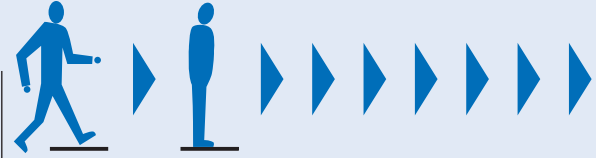
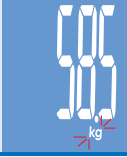
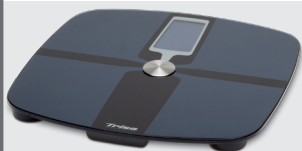
DE | FR | IT | EN | ES

**Table 1:** Guidance and MANUFACTURER'S declaration – ELECTROMAGNETIC EMISSIONS for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacturer's declaration – electromagnetic emissions		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 2	The device must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

**Table 2:** Guidance and MANUFACTURER'S declaration – electromagnetic IMMUNITY – for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacturer's declaration – electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	$\pm 6$ kV contact $\pm 8$ kV air	$\pm 6$ kV contact $\pm 8$ kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	$\pm 2$ kV for power supply lines $\pm 1$ kV for input/output lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	$\pm 1$ kV line(s) to line(s) $\pm 2$ kV line(s) to earth	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	$<5\% U_T$ ( $>95\%$ dip in $U_T$ ) for 0.5 cycle	Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery.
	$40\% U_T$ (60% dip in $U_T$ ) for 5 cycles		
	$70\% U_T$ (30% dip in $U_T$ ) for 25 cycles		
	$<5\% U_T$ ( $>95\%$ dip in $U_T$ ) for 5 s		
Power frequency (50/60Hz) magnetic field IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE $U_T$ is the a.c. mains voltage prior to application of the test level.			

**Table 4:** Guidance and MANUFACTURER'S declaration – electromagnetic IMMUNITY – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacturer's declaration – electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
			<p>Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p><b>Recommended separation distance</b></p>
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	Not applicable	$d = \left[ \frac{3.5}{V_1} \right] \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	80 MHz to 800 MHz: $d = 1.167 \sqrt{P}$ 800 MHz to 2.5 GHz: $d = 2.333 \sqrt{P}$
			<p>where <math>P</math> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <math>d</math> is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey<sup>a</sup>, should be less than the compliance level in each frequency range<sup>b</sup>.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<b>Note 1</b>	At 80 MHz and 800 MHz, the higher frequency range applies.		
<b>Note 2</b>	These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		
<sup>a</sup>	Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.		
<sup>b</sup>	Over the frequency range 150 kHz to 80 MHz, field strengths should be less than $[V_1]$ V/m.		

**Table 6:** Recommended separation distances between portable and mobile RF communications equipment and the ME EQUIPMENT or ME SYSTEM – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Recommended separation distances between portable and mobile RF communications equipment and the device.			
The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = \left\lceil \frac{3,5}{V_1} \right\rceil$	80 MHz to 800 MHz $d = 1.167$	800 MHz to 2.5 GHz $d = 2.2333$
$\sqrt{P}$			
0.01	Not applicable	0.117	0.233
0.1	Not applicable	0.369	0.738
1	Not applicable	1.167	0.233
10	Not applicable	3.690	7.338
100	Not applicable	11.67	23.33
For transmitters rated at a maximum output power not listed above, the recommended separation distance <b>d</b> in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where <b>P</b> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
<b>Note 1</b>	At 80MHz and 800MHz, the separation distance for the higher frequency range applies.		
<b>Note 2</b>	These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.		

## Garantie-Hinweis | Informations de garantie | Dichiriazione di garanzia | Warranty information | Garantía – Nota

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**EN** With this appliance you get a 2 year guarantee starting from the date of purchase. Under the guarantee, in the case of material or manufacturing defects, the appliance will be replaced or repaired. Replacement with a new appliance or money back return is not possible. Excluded from the guarantee is normal wear and tare, commercial use, alterations to the appliance as purchased, cleaning activities, consequences of improper use or damage by the purchaser or a third person, damage, that can be attributed to external circumstances or caused by the batteries. The guarantee requires that the faulty appliance be returned at the purchaser's expense along with the retail outlets dated and signed guarantee card or a sales receipt.

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**Switzerland** Trisa Electronics AG
Satellitentechnik Vertriebs GmbH
HKS-GmbH@T-online.de
Service GmbH