# beurer

## MG 850

fascia releaZer

|           |   |     |
|-----------|---|-----|
| <b>DE</b> | Faszien-Massagegerät & Deep-ReleaZer Aufsätze<br>Gebrauchsanweisung .....                     | 2   |
| <b>EN</b> | Fascia massager & Deep-ReleaZer attachments<br>Instructions for use .....                     | 23  |
| <b>FR</b> | Appareil de massage du fascia et embout Deep-ReleaZer<br>Mode d'emploi .....                  | 43  |
| <b>ES</b> | Aparato de masaje para fascias y accesorios Deep-ReleaZer<br>Manual de instrucciones .....    | 65  |
| <b>IT</b> | Apparecchio per il massaggio fasciale e accessori Deep-ReleaZer<br>Istruzioni per l'uso ..... | 86  |
| <b>TR</b> | Bağdoku masaj cihazı ve Deep-ReleaZer başlıkları<br>Kullanım kılavuzu .....                   | 107 |
| <b>RU</b> | Прибор для массажа фасций с насадками Deep-ReleaZer<br>Инструкция по применению .....         | 128 |
| <b>PL</b> | Aparat do masażu powięziowego i nasadki Deep-ReleaZer<br>Instrukcja obsługi .....             | 150 |

# DEUTSCH

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# ENGLISH

**Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.**

## **Warning**

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to adhere to this instruction will void the warranty.

## Dear customer,

thank you for choosing a product from our range. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty and air. Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

With kind regards,  
Your Beurer team

## Table of contents

|   |           |  |
|---|-----------|--|
| <b>1. Included in delivery .....</b>  | <b>25</b> |  |
| <b>2. Signs and symbols .....</b>   | <b>25</b> |  |
| <b>3. Warnings and safety notes .....</b>                                   | <b>27</b> |  |
| <b>4. Getting to know your device ....</b>                                  | <b>28</b> |  |
| <b>5. Intended use .....</b>  | <b>30</b> |  |
| <b>6. Unit description .....</b>  | <b>32</b> |  |
| 6.1 Unit description for fascia<br>massager .....                           | 32        |  |
| 6.2 Unit description for<br>Deep-ReleaZer attachments..                     | 33        |  |
| <b>7. Initial use .....</b>   | <b>34</b> |  |
| <b>8. Use .....</b>   | <b>35</b> |  |
| 8.1 Preparing for the application   | 35        |  |
| 8.2 Switching on the Fascia-ReleaZer®<br>/ selecting the vibration setting. | 36        |  |
| 8.3 Mounting the attachments....  | 37        |  |
| 8.4 Exercises .....   | 38        |  |
| <b>9. Cleaning and maintenance .....</b>                                    | <b>39</b> |  |
| <b>10. Accessories and replacement<br/>    parts.....</b>                   | <b>39</b> |  |
| <b>11. Disposal .....</b>   | <b>39</b> |  |
| <b>12. Technical specifications .....</b>                                   | <b>40</b> |  |
| <b>13. Warranty/service .....</b>   | <b>41</b> |  |

## 1. Included in delivery

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

- 1 x Fascia-ReleaZer®
- 1 x set of Deep-ReleaZer attachments
- 1 x micro USB charging cable
- 1 x mains adapter
- 1 x training DVD with training exercises
- 1 x exercise overview

## 2. Signs and symbols

The following symbols appear in these instructions for use:

|   |                     |   |
|---|---------------------|---|
|  | <b>WARNING</b>      | Warning notice indicating a risk of injury or damage to health. |
|  | <b>IMPORTANT</b>    | Safety note indicating possible damage to the device/accessory. |
|  | <b>Note</b>         | <b>Product information</b><br>Note on important information     |
|  | For indoor use only |   |

|   |  |
|---|--|
|   | <p>Read the instructions</p>   |
|   | <p>Dispose of packaging in an environmentally friendly manner</p>  |
|  | <p><b>Protection class II device</b><br/>The device is double-protected and thus corresponds to protection class 2</p>     |
|  | <p>Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE</p>                       |
|  | <p>Do not dispose of batteries containing hazardous substances with household waste.</p>                                   |
|  | <p>Manufacturer</p>  |
|  | <p><b>CE labelling</b><br/>This product satisfies the requirements of the applicable European and national directives.</p> |
|  | <p>The products demonstrably meet the requirements of the Technical Regulations of the EAEU.</p>                           |

|   |   |
|---|---|
|   | Power supply unit, switch mode; switch mode power supply unit; SMPS |
|  | Safety isolating transformer, short-circuit proof                   |
|   | Polarity  |
|  | Direct current  |
|  | Energy efficiency class 6   |

### 3. Warnings and safety notes

#### WARNING

#### Safety note regarding use of the Fascia-ReleaZer®

- Keep the packaging material away from children. Choking hazard!
- The Fascia-ReleaZer® and the attachments are not medical devices. They are massage treatment devices.
- Before use, ensure that the device and accessories show no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.
- Repairs must only be carried out by Customer Services or authorised suppliers. Under no circumstances should you open or repair the Fascia-ReleaZer® yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to comply will void the warranty. Never attempt to repair the device yourself!

- Do NOT use the device on the head.
- Do NOT use the device on the face.
- Do NOT use the device on the genitals.
- Only use the device with the supplied accessories.

## **Notes on handling batteries**

- If your skin or eyes come into contact with fluid from the battery cell, flush out the affected areas with plenty of water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush the batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter Initial use).
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

## **4. Getting to know your device**

### **What are fascia?**

Fasciae are fine, tough connective tissue fibres that pass through our entire body like a three-dimensional network. Human fasciae surround all muscles, bones, organs and nerve cells in the body and ensure that when we move, the individual parts of our body glide past one another and remain in place.

Fasciae also act as links between the muscle fibres and ensure that we are able to use our muscle strength. The strain energy between the tendons and bones is stored and then suddenly released upon relaxing.

Fasciae shorten and harden if they are not used on a regular basis (similarly to muscles). The structures of the fasciae stick together and thicken, which causes them

to become porous and lose their elasticity. This results in small tears and inflammation, which also affects the transmission of force and your coordination. In isolated cases, thickened fasciae may also press down on the nerves and cause discomfort. There may be a variety of reasons as to why fasciae have shortened or stuck together, such as a lack of movement, protective postures or psychological stress.

## **How does the Fascia-ReleaZer® work?**

If fascia tissue has become stuck together or hardened, you can improve its suppleness and elasticity by alternately applying pressure and tension. By using the Fascia-ReleaZer® as a self-help device, you can positively stimulate these receptors in our tissue. This will increase and speed up the regenerative capacity of the tissue. Your blood supply will also be improved and your nerves will function better. The stimulation from the Fascia-ReleaZer® allows new nutrients to reach the cells and makes it easier to transport away waste. In this way, the tissue is replaced and the elasticity and resilience of the tissue in the treated area are improved.

### **According to study:**

Randomised, Double Controlled, Standardised, Clinical Study, Europe.

### **Performed by:**

Center for Integrative Therapy, Germany.

Institute of Medical Psychology and Behavioral Neurobiology, University of Tübingen, Germany. Research Institute on Health Sciences (IUNICS), University of Balearic Islands, Palma, Spain.

## **How do the Deep-ReleaZer attachments work?**


The various attachments help to relieve tension in the back region. They also relax the muscles in the back and neck area, as well as the hips and thighs. The Fascia-ReleaZer® can be used with or without attachments.

Another way to boost the elasticity of the tissue is with the vibration massage, which relaxes the muscles. This reduces stiffness, which in turn improves your movement.

The massage treatment must feel pleasant at all times. A slight haematoma may well occur – this is a positive therapeutic effect as the use of the Fascia-ReleaZer® can selectively encourage your blood flow. In isolated cases areas of redness may also occur on the skin. If you experience excessive skin irritation, stop using the Fascia-ReleaZer® and consult a doctor.

## When are the different attachments used?

The "soft" attachment is used to start the application and for slight muscle pain.

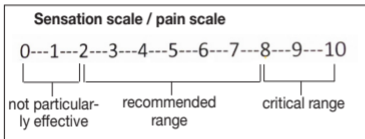After a few uses, you can move on to the "medium" attachment. This strengthens the muscles. At the next level, the "hard" attachment can be used to intensify application.

To ensure that the treatment is pleasant and tailored to you, we

recommend the use of the sensation/pain scale. This will help you to find the ideal level of pressure for you. Since each person has a different perception of pain, we recommend a massage treatment in the recommended range (2-8). If the sensation rating lies between 0 and 2, the treatment will not be particularly effective as the Fascia-ReleaZer® incl. attachments is unable to penetrate deep into the tissue. We would advise against a massage treatment in the critical range (more than 8), since excessive pressure is exerted on the tissue in this range.

## 5. Intended use

### WARNING

Only use the Fascia-ReleaZer® and the attachments for massaging human muscle fasciae, backs, necks, hips and arms.

The Fascia-ReleaZer® and the attachments are not medical devices, but massage treatment devices.

The device must only be used for the purpose for which it is designed and in the manner specified in the instructions for use. Improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or careless use.

Use the Fascia-ReleaZer® and the attachments

- Not with changes caused by disease or injury in the back region (e.g. disc problems, open wounds),
- not during pregnancy,
- never while you sleep,
- not in the car,
- not on animals,
- not during activities during which an unexpected reaction may be dangerous,
- not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol),
- never longer than 15 minutes (risk of overheating) and allow it to cool off for at least 15 minutes before using it again.


Please consult your doctor if you are unsure whether or not the Fascia-ReleaZer® is suitable for you.

You should consult your doctor before using the Fascia-ReleaZer® if:

- you suffer from a serious illness or are recovering from an operation,
- you suffer from any of the following: common cold accompanied by a fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica) or acute inflammations.
- you are fitted with a cardiac pacemaker, implant or other aid,
- in case of pains of unclarified origin.

## 6. Unit description


### 6.1 Unit description for fascia massager
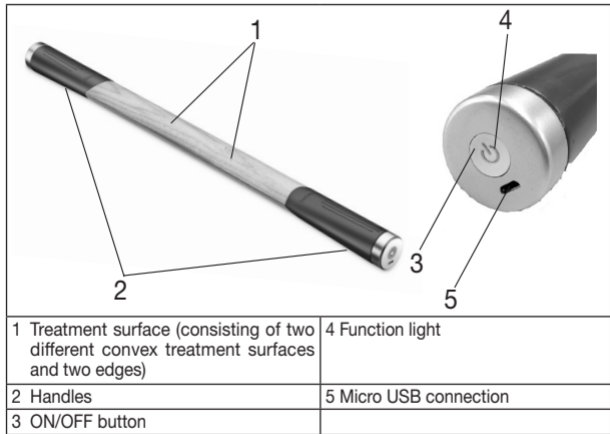
| Function light status  | Vibration function      |
|--|-------------------------|
| Function light shows steady green  | = Static vibration      |
| Function light flashes green   | = Wave-action vibration |

| Function light status  | Charging function (charging cable plugged in) |
|---|---|
| Function light flashes yellow   | = Charging                                    |
| Function light shows steady yellow  | = Charging complete                           |
| Function light flashes red  | = Low battery                                 |

## 6.2 Unit description for Deep-ReleaZer attachments
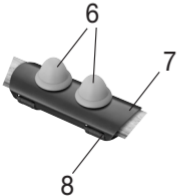
|  |   |
|--|---|
|  | 6 Attachments (3 degrees of hardness: soft/medium/hard) |
|  | 7 Holder for attachments                                |
|  | 8 Hook-and-loop fastener for securing the holder        |
|  |   |


## 7. Initial use

### Note

Before using the Fascia-ReleaZer® for the first time, let it charge for a minimum of 2.5 hours. Proceed as follows:

- Connect the charging cable to the mains adapter and the Fascia-ReleaZer® (*Figure 1*). Then insert the mains adapter into a suitable socket. Alternatively, you can also charge the Fascia-ReleaZer® with the charging cable via the USB port on your computer or laptop (*Figure 2*).

Whilst charging, the function light flashes yellow. When the Fascia-ReleaZer® is fully charged, the function light shows steady yellow. The battery will last for approx. 11 hours of operation when fully charged. Charge the Fascia-ReleaZer® again as soon as the function light flashes red during operation.
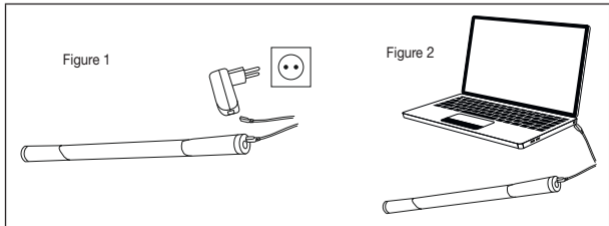
## 8. Use

### 8.1 Preparing for the application

We recommend treatment with or without plant-based massage oil. If you have thick body hair, we recommend treatment with plant-based massage oil. However, we generally advise using the device on smooth, hair-free skin. Do not use any oils or creams with a strong cooling or heating effect, as they may cause skin irritation.

#### **Treatment without massage oil: (recommended if you have no or little body hair)**

A massage treatment without massage oil provides a surface massage. This massage treatment activates the reflex zones and treats the surface of the tissue.

#### **Treatment with massage oil: (recommended if you have thick body hair or if you would like a deep massage)**

If you opt for massage oils, the massage will penetrate deep into the tissue. This will strongly promote your circulation.

#### **IMPORTANT**

Only use plant-based oils as other oils may stain or damage the wood surface of the Fascia-ReleaZer®.

1. Put a couple of drops of plant-based oil in the palm of your hand.
2. After the massage treatment, wipe the wood surface of the Fascia-ReleaZer®, the attachments and the treated area of the body with a towel
3. After the massage treatment, wipe the wood surface on the Fascia-ReleaZer® and the treated area with a cloth.

The amount of oil to be used depends on the effect you wish to achieve. More oil causes a faster movement and therefore a more intensive reflex action and circulation of blood to the skin.

Less oil causes a slower treatment that penetrates deeper, with a very relaxing effect that loosens the tissue.

It is also possible to combine the two. To do so, apply a larger amount of oil onto the area to be treated and move the Fascia-ReleaZer® over the area quickly. Then wipe off some of the oil with a cloth and treat the same area again, but this time with slow movements.

## **8.2 Switching on the Fascia-ReleaZer® / selecting the vibration setting**

Press the ON/OFF button to switch on the Fascia-ReleaZer®. The function light shows steady green and the first vibration setting is activated.

| <b>Action</b>           |   | <b>Vibration setting</b>  |
|-------------------------|---|---|
| 1st press of the button | = | The Fascia-ReleaZer® switches on and the function light shows steady green.<br>Static vibration |
| 2nd press of the button | = | The function light flashes green.<br>Wave-action vibration                                      |
| 3rd press of the button | = | The Fascia-ReleaZer® switches off and the function light goes out.                              |

The Fascia-ReleaZer® can be used statically (targeted application) or dynamically (moving application).

### **Static application:**

This method loosens specific areas of hardness. Press down evenly on the same spot approx. 3 times in succession.

### **Dynamic application:**

The movement smooths out the fasciae. Moving the device quickly forwards and backwards provides a surface massage for the upper tissue. Moving the device more slowly allows the massage effect to penetrate deeper.

Avoid striking the device against the treatment area as this may result in haematoma.

### **Treatment surfaces:**

We recommend starting the massage treatment with one of the two different convex treatment surfaces. After approx. one week, the massage treatment can be continued with both edges. The edges penetrate deep into the tissue and provide a deep massage treatment.

## **8.3 Mounting the attachments**

Secure the attachments you wish to use (soft/medium/hard) to the holder. To do so, twist them in a clockwise direction onto the holder. Then open the hook-and-loop fastener and pull the holder over the fasciae massager. The attachment should be placed in the middle of the device. Then close the hook-and-loop fastener.

## 8.4 Exercises

### **Note**

We recommend treating each area (lower body, mid-body and upper body) three times a week for 15 minutes each time.

When carrying out a session, we recommend training all three areas. However, you may also choose to focus on one of the above areas.

After each massage treatment, wipe down the Fascia-ReleaZer® and the attachments with a dry cloth or hand towel. Ensure that no traces of oil that could attack the surface are left on the device.

**All exercises can be found on our enclosed training DVD and in our exercise overview.**

### **Fascia massager**

#### **Exercises for the lower body**

Lower leg, thighs, shin, calves

#### **Exercises for the mid-body**

Lower back (not on the spine), bottom

#### **Exercises for the upper body**

Neck, upper and lower arms

#### **Deep-ReleaZer attachments**

Back, neck, hip, thigh, inner arm and palm

## 9. Cleaning and maintenance

The service life of the Fascia-ReleaZer® depends on careful use.

### IMPORTANT

- Protect the Fascia-ReleaZer® and the attachments from knocks, dust, chemicals, drastic changes in temperature, electromagnetic fields and nearby sources of heat (ovens, heaters).
- Clean the Fascia-ReleaZer® and the attachments with a dry cloth after every use. Do not use any abrasive cleaning products.
- Cleaning must not be performed by children unless supervised.

## 10. Accessories and replacement parts

Replacement parts and wearing parts are available from the corresponding service address (according to the service address list). Please state the corresponding order number.

|                             | Item number and/or order number |
|-----------------------------|---------------------------------|
| Attachments replacement set | 163.819                         |

## 11. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life.

Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

## Rechargeable battery disposal

The empty, completely flat rechargeable batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the rechargeable batteries.

The codes below are printed on rechargeable batteries containing harmful substances:

Pb = Battery contains lead,

Cd = Battery contains cadmium,

Hg = Battery contains mercury.

## 12. Technical specifications

|                                     |  |
|-------------------------------------|--|
| Mains operation:<br>Input<br>Output | 100–240 V~, 50/60 Hz, 0.5 A<br>5.0 V <u>===</u> 1.0 A, 5.0 W |
| Average active efficiency:          | ≥ 77.11 %  |
| No-load power consumption:          | ≤ 0.065 W  |
| Dimensions                          | 60 x 4.7 x 3.7 cm  |
| Weight: approx                      | approx. 600g   |

|  |                             |
|--|-----------------------------|
| <b>Battery:</b><br>Capacity<br>nominal voltage<br>type designation | 2200 mAh<br>3.7 V<br>Li-Ion |
|--|-----------------------------|

### 13. Warranty/service

Beurer GmbH, Söflinger Straße 218, 89077 Ulm, Germany (hereinafter referred to as “Beurer”) provides a warranty for this product, subject to the requirements below and to the extent described as follows.

**The warranty conditions below shall not affect the seller’s statutory warranty obligations which ensue from the sales agreement with the buyer. The warranty shall apply without prejudice to any mandatory statutory provisions on liability.**

Beurer guarantees the perfect functionality and completeness of this product.

The worldwide warranty period is 3 years, commencing from the purchase of the new, unused product from the seller.

The warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use. German law shall apply.

During the warranty period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, Beurer shall carry out a repair or a replacement delivery free of charge, in accordance with these warranty conditions.

**If the buyer wishes to make a warranty claim, they should approach their**

**local retailer in the first instance: see the attached “International Service” list of service addresses.**

The buyer will then receive further information about the processing of the warranty claim, e.g. where they can send the product and what documentation is required.

A warranty claim shall only be considered if the buyer can provide Beurer, or an authorised Beurer partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this warranty:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by Beurer;
- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the warranty period under any circumstances.

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Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany