
beurer

SC 30


DE	Elektrische Zahnbürste	
	Gebrauchsanweisung	6
EN	Electric toothbrush	
	Instructions for use.....	13
FR	Brosse à dents électrique	
	Mode d'emploi.....	19
ES	Cepillo de dientes eléctrico	
	Manual de instrucciones	25
IT	Spazzolino elettrico	
	Istruzioni per l'uso	32
TR	Elektrikli diş fırçası	
	Kullanım kılavuzu	38
RU	Электрическая зубная щетка	
	Инструкция по применению	44
PL	Elektryczna szczoteczka do zębów	
	Instrukcja obsługi	51
NL	Elektrische tandenborstel	
	Gebruiksaanwijzing	58
DA	El-tandbørste	
	Betjeningsvejledning.....	64
SV	Eltandborste	
	Bruksanvisning.....	70
NO	Elektrisk tannbørste	
	Bruksanvisning.....	76
FI	Sähköhammasharja	
	Käyttöohje.....	82

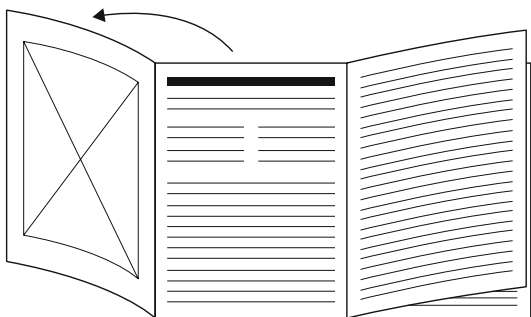
- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Unfold page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- ES** Despliegue la página 3 antes de leer las instrucciones de uso.
- IT** Prima di leggere le istruzioni per l'uso aprire la pagina 3.
- TR** Kullanım kılavuzunu okumadan önce 3. sayfayı açın.
- RU** Перед чтением инструкции по применению разложите страницу 3.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit om de gebruiksaanwijzing te kunnen lezen.
- DA** Fold side 3 ud, før du læser betjeningsvejledningen.
- SV** Vik ut sidan 3 innan du läser bruksanvisningen.
- NO** Åpne side 3 før du leser bruksanvisningen.
- FI** Käännä sivu 3 auki ennen käyttöohjeen lukemista.

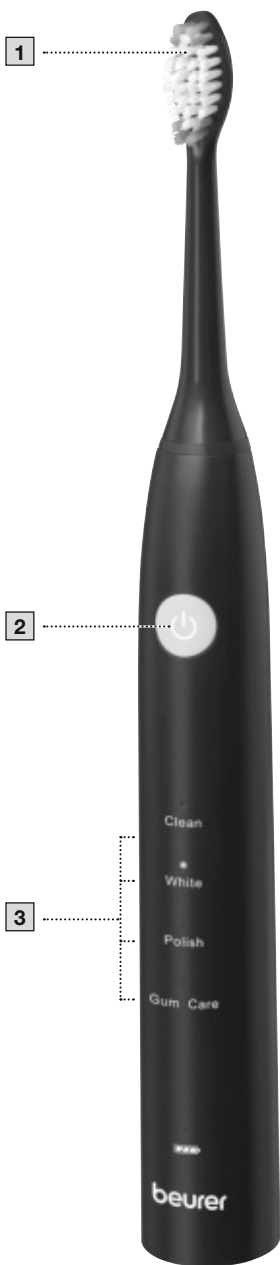
A


B

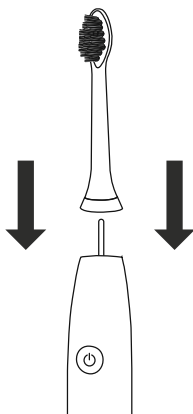
C

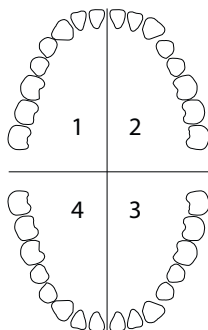D

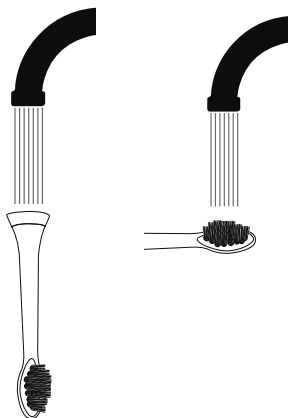E

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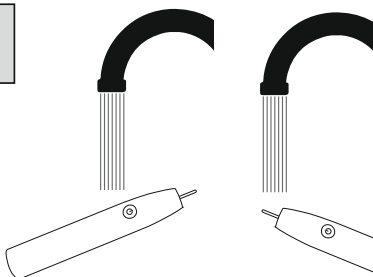Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the next user with these instructions for use.

CONTENTS

1. Warnings and safety notes	13
2. Included in delivery	14
3. Signs and symbols	15
4. Intended purpose	16
5. Device description.....	16
6. Usage	16
6.1 Preparing for use.....	16
6.2 Starting use.....	17
6.3 Travel lock	17
7. Cleaning and maintenance.....	17
8. Accessories and replacement parts.....	17
9. Disposal.....	18
10. Technical specifications	18
11. Guarantee.....	18

1. WARNINGS AND SAFETY NOTES

⚠ WARNING

- Keep the device and the charging cable away from water or other liquids.
- Do not place the device near a sink or bathtub full of water.
- This device may be used by children over the age of 8 years and by people with reduced physical, sensory or mental capabilities or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the risks of using the device.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.

- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- If the mains connection cable is broken, you must dispose of the charging cable.
- Only use the USB-C charging cable supplied with your device.
- Never use the device with the mains adapter near a bath, sink, shower or other containers filled with water or other liquids – risk of electric shock.
- The device may only be charged with a mains adapter (not included in delivery) with a safety extra-low voltage output of 5 V and max. 1 A.
- Under no circumstances should you open or repair the device yourself, otherwise proper functionality can no longer be guaranteed. Failure to comply will invalidate the warranty.

Notes on handling batteries

⚠ WARNING

- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- Risk of explosion! Do not throw batteries into a fire.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Do not disassemble, split or crush the batteries.

⚠ CAUTION

- Protect batteries from excessive heat.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter “Preparing for use”).
- In order to achieve as long a battery service life as possible, fully charge the battery at least 2 times a year.

2. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been

removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 electric toothbrush
- 2 brush heads
- 1 USB-C charging cable
- Instructions for use

3. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	DANGER
Indicates an imminent danger. If it is not avoided, death or serious injury will occur.	
	WARNING
Indicates a potentially imminent hazard. If it is not avoided, there is a risk of death or serious injury.	
	CAUTION
Indicates a potentially imminent hazard. If it is not avoided, slight or minor injuries may result.	
	NOTICE
Indicates a potentially harmful situation. If it is not avoided, the device or something in its vicinity may be damaged.	
	Product information Note on important information
	Read the instructions
	Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).
	Do not dispose of batteries containing hazardous substances with household waste
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	Standard for independent certification of recycled material
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = Material code, B = Material number: 1–7 = Plastics, 20–22 = Paper and cardboard
	Conformity assessed marking for the United Kingdom
	Protection class III device

	Importer symbol
	Manufacturer
	Protection against temporary submersion. Device can be immersed to a depth of 1 metre for up to 30 minutes.
	USB / Micro-USB connection
	For indoor use only
	Direct current
	Packaging made from recycled material

4. INTENDED PURPOSE

Only use the device to clean teeth on humans. Do not use the device on animals! The device is not suitable for children under 8 years. The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

- 1** Brush head
- 2** ON/OFF/Mode button
- 3** LED display
 - 4 modes: Clean, White, Polish, Gum Care
 - Battery level
- 4** USB-C charging cable

6. USAGE

6.1 Preparing for use

Charging the electric toothbrush **A**

Fully charge the electric toothbrush before first use. To do so, proceed as follows:

1. Connect the USB-C charging cable to the electric toothbrush. Plug the USB-C charging cable into a USB port or into a suitable mains adapter and plug it into a socket.
2. Whilst charging, the LED pulsates red. As soon as the electric toothbrush is fully charged, the LED lights up green continuously.

Fitting the brush head **B**

Place the brush head **1** onto the device from above until it clicks into place. Make sure that the brush head **1** points in the same direction as the ON/OFF/Mode button **2**.

Indicator bristles

Indicator bristles lose their colour over time depending on how each person uses the device, thus indicating the best time to change the brush head **1**.

We recommend changing the brush head **1** after 3 months. If the bristles of the brush head are worn or frayed, they may need to be replaced more frequently.

6.2 Starting use

1. To switch on the electric toothbrush, briefly press the ON/OFF/Mode button **2**. The electric toothbrush pauses briefly every 30 seconds to remind you to clean all 4 jaw quadrants for the same amount of time **C**. The electric toothbrush has an automatic timer and vibrates two times after 2 minutes. The toothbrush switches off automatically after 6 minutes.
2. Select your desired mode by briefly pressing the ON/OFF/Mode button **2**. The currently selected mode is shown accordingly by the corresponding mode LED. The following modes are available:

Clean	For daily cleaning
White	For whiter teeth
Polish	For removing superficial stains
Gum Care	For looking after the gums

3. To cancel the application prematurely, briefly press the ON/OFF/Mode button **2**.

6.3 Travel lock

The travel lock prevents the device from being switched on by accident (e.g. on trips, when being carried in a suitcase).

1. To activate the travel lock, press and hold the ON/OFF/Mode button **2** for 5 seconds. During activation, the LED flashes green 5 times.
2. If you press the ON/OFF/Mode button **2** briefly when the device is locked, the LED will also flash green.
3. To deactivate the travel lock, press and hold the ON/OFF/Mode button **2** again for 5 seconds. During deactivation, the LED lights up green continuously. The LED goes out as soon as the travel lock has been deactivated. The device starts up.

7. CLEANING AND MAINTENANCE

Clean the device and the brush head **1** after each use. To do so, proceed as follows:

1. Pull the brush head **1** upwards and off.
2. Clean the brush head **1** under running water. Make sure to clean the inside of the brush head **1** as well. **D**
3. Clean the device under running water. However, do not submerge the device in water. Water might enter the device and damage it. **E**
4. Regularly clean the charger **5** with a damp cloth and a mild soap solution. Do not use any abrasive cleaning products. Keep the inductive charger **5** and the USB charging cable with mains adapter away from water.

8. ACCESSORIES AND REPLACEMENT PARTS

To purchase accessories and replacement parts, please visit www.beurer.com or contact the corresponding service address (according to the service address list) in your country. Accessories and replacement parts are also available from retailers.

Designation	Item number and/or order number
Classic Clean brush heads 4-pack	105.21
Premium Whiten brush heads 4-pack	105.31
Premium Plaque Eraser brush heads 4-pack	105.22
Premium Protection Plus brush heads 4-pack	105.23

9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local authorities.

The empty, completely flat batteries must be disposed of in specially designated collection boxes, at recycling points or at electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

Removing the battery

- 1. Use a screwdriver to open the electric toothbrush. Insert the screwdriver into the USB-C port.
- 2. Now apply a little pressure and turn anti-clockwise. The end cap on the toothbrush can now be removed.
- 3. Fold the bottom of the toothbrush to the side.
- 4. Pull out the inner housing of the toothbrush.
- 5. Remove the battery from the inner housing.
- 6. Dispose of the battery separately in specially designated collection boxes, at recycling points or at electronics retailers.

10. TECHNICAL SPECIFICATIONS

Type	SC 30
Dimensions	Device: 183 x 29 mm
Weight	Device: 119 g (toothbrush + 2 brush heads)
Battery: Capacity Nominal voltage Type designation	600 mAh 3.7 V Lithium-ion battery
Power supply	5.0 V $\overline{\text{---}}$ 1.0 A, 5.0 W

Subject to technical changes

11. GUARANTEE

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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Beurer GmbH • Söflinger Str. 218
BEURER UK Ltd., Suite 16,