SKG


Waist Massager
VS500

USER MANUAL

Catalog of Instruction Manuals

English User Manual.....	1
Manuel de l'utilisateur en français.....	7
Deutsches Benutzerhandbuch.....	13
Manuale d'uso italiano.....	19
Manual del usuario en español.....	25
Polska instrukcja obsługi.....	31
Посібник користувача українською мовою.....	37
Türk Kullanıcı Elçisi.....	43

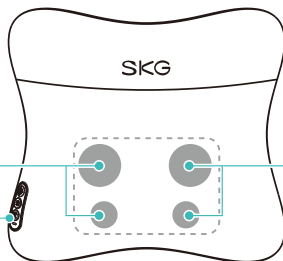Massage Head
/Hot Compress

Control
Panel

Massage Head
/Hot Compress


On/Off button

(Press and hold to turn on/off;
Short press to switch mode)


Hot compress button on the front
(Short press to turn on/off hot
compress)

Hot compress button on the back
(Short press to switch intensity level)

Charging port

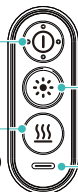Hot compress area

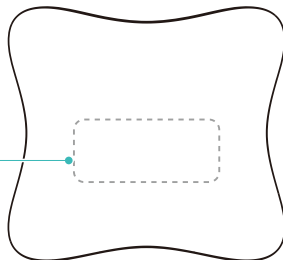(The diagram above is for function illustration purposes only and may differ in appearance from the actual product)

Note: Please do not massage the same area for too long.

1. This product features a massage and hot compress area on one side, and a hot compress area on the other for leaning on.

2. All operation buttons are located on the control panel of the massager. To power on or off, press and hold the power button on the control panel for 3 seconds until the massage head starts working; to switch between massage modes, press the power button briefly; to turn on or off the hot compress function, press the Hot Compress button on the front briefly; to activate the hot compress function on the back, press the back Hot Compress button briefly, and press it again to cycle through intensity levels or press it for the third time to turn off the hot compress function.

3. A solid green power indicator light signifies that the battery is fully charged, while an orange flashing light indicates that the battery is about to run out. When the product is out of power, it can be charged using the Type-C charging port. The light on the charging port will indicate different charging statuses: an orange breathing light indicates that the device is charging, and a green breathing light indicates that the battery is fully charged.

4. To end the massage, press and hold the power button on the control panel for 3 seconds to turn off the massage function.

Basic Parameters

Product name	Waist Massager		
Model	VS500	Rated input voltage	5V $\equiv$
Battery capacity	3200mAh/23.68Wh	Rated input power	10W

This product is rechargeable, and users should charge it with a Type-C charging cable. Chargers used must have an output not more than 15W.

- Please never use the device with the following medical electronic devices: implantable medical electronic devices such as cardiac pacemakers, life-sustaining medical electronic devices such as artificial heart-lung machine, medical electronic devices such as electrocardiograph
- This product is a personal care massager product and is not intended to replace pharmaceuticals and medical devices.
- Pregnant and postpartum women should use this product with caution and only under the guidance and permission of a physician.
- Individuals with sensitive skin should use this product only with the approval and guidance of a physician.
- This product contains a lithium battery. Do not disassemble or modify the product, immerse it in water, or throw it in fire.
- When using the product, the massage head should not contact with the following areas: waist wounds and scars, areas with internal malignant tumors, areas with internal metal implants, and areas of acute purulent inflammation.
- The following individuals are prohibited from using this product: those with coagulation disorders, those with spastic paralysis, those in the spastic phase following a cerebrovascular accident, those with severe cognitive impairments, those with nervous irritability, those with contractures or deformities of joints, infants and children who have lost all or part of their consciousness.
- The following individuals should consult a physician before using the product: patients who are under medical treatment and feel unwell, patients with infectious diseases, patients with a body temperature exceeding 38 degrees Celsius (feverish), patients with abnormal blood pressure, patients with cardiac or neurological abnormalities, patients with sensory impairments, patients who cannot freely express their own will, patients with osteoporosis, and patients with spinal fractures.
- In case of any discomfort during use, please stop the product immediately.
- The battery level may deplete quickly at the initial use or after the product has not been used for an extended period. Please charge the battery following the User Manual to ensure its normal operation.
- Do not plug or unplug the charger with wet hands, as this may cause electric shocks or accidents.
- Recommended for use by people of 16 years and older.
- If the massager is not used for a long time, it is recommended to charge it once a month to keep the battery level above 20%.

Please support the recycling of rechargeable batteries (lithium-ion batteries)

 Danger	The battery used in this product is a special rechargeable battery, which should not be used in other products. Do not charge the rechargeable battery that has been taken out. • Do not overheat the battery by throwing it into the fire. • Do not pierce with nails, impact by gravity or disassemble and modify. • Do not connect the battery cathode and anode directly with a metal. • Do not store it together with necklaces, hairpins, etc. • Do not charge, use and place the battery in a hot place, such as by the fireplace or under the hot sun The above situations can all lead to heat, ignition, and rupture.
 Warning	Please store the removed rechargeable battery out of reach of infants and young children to prevent accidental ingestion. • If ingestion occurs, seek immediate medical attention. • In the event of battery leakage, do not touch the liquid with bare hands and follow these steps for handling. • When the battery fluid gets in your eyes, there is a risk of blindness. Do not rub your eyes directly with your hands, but flush them with clean water immediately and seek medical treatment. When the battery fluid is in contact with skin or clothes, it may cause skin inflammation or injury. Please flush with clean water, and seek medical treatment immediately if you feel uncomfortable. • Disassembly of the product is only permitted when it is disposed of. Failure to comply may result in accidents such as fires, electric shocks, or injuries. If repair is needed (excluding battery replacement), please consult the store where you purchased the product.

The lithium-ion battery used in this product is a valuable recyclable resource. Therefore, please hand it over to a professional organization that can conduct recycling and regeneration of rechargeable batteries when it reaches the end of its life cycle and will be discarded. The battery must be removed from the massager before discarding. The massager must be powered off before the battery is removed. The battery should be disposed of safely.

Steps to remove the battery:

1. Ensure the product is powered off.
2. Open the zipper of the fabric cover.
3. Loosen the screws on the gear case cover.
4. Disconnect the battery output wire terminals.
5. Remove the unit.
6. Remove the battery.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device in compliance with the essential requirements and other relevant provisions of Directive (EU)2023/988.

CAUTION : RISK OF EXPLOSION IF BATTERY IS REPLACED BY AN INCORRECT TYPE. DISPOSE OF USED BATTERIES ACCORDING TO THE INSTRUCTIONS.

EU Compliance Statement:

SKG Health technology Co.,Ltd. Here by declares that this device is in compliance with the essential requirements and other relevant provisions of the Directive (EU)2023/988. A copy of the EU Declaration of Conformity is available online at <https://www.skg.com/pages/skg-declaration-of-conformity>.

Old electrical appliances must not be disposed of together with the residual waste, but have to be disposed of separately. The disposal at the communal collecting point via private persons is for free. The owner of old appliances is responsible to bring the appliances to these collecting points or to similar collection points. With this little personal effort, you contribute to recycle valuable raw materials and the treatment of toxic substances.

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AB uygulama ifadesi:
SKG Health Technologies Co., Ltd.