



KERN & Sohn GmbH

Ziegelei 1
72336 Balingen-Frommern
Germany

+0049-[0]7433-9933-0
+0049-[0]7433-9933-149
info@kern-sohn.com

Installation Instructions Stands

KERN SXC, IXC

Type

TSXC-A01-A
TIXC-A01-A
TIXC-A02-A

Version 1.0
2023-06
GB



TSXC-A01-A_TIXC-A01-A_TIXC-A02-A-IA-e-2310



Stands
Version 1.0 2023-06
Installation Instructions

Contents


1 Technical data..... 3
1.1 Dimensions 3
2 Scope of delivery 6
3 Installation..... 7

1 Technical data




KERN	TSXC-A01-A	TIXC-A01-A	TIXC-A02-A
Stand height (mm)	330	330	600

1.1 Dimensions




- Dimensions in mm
- Drawings on the following page

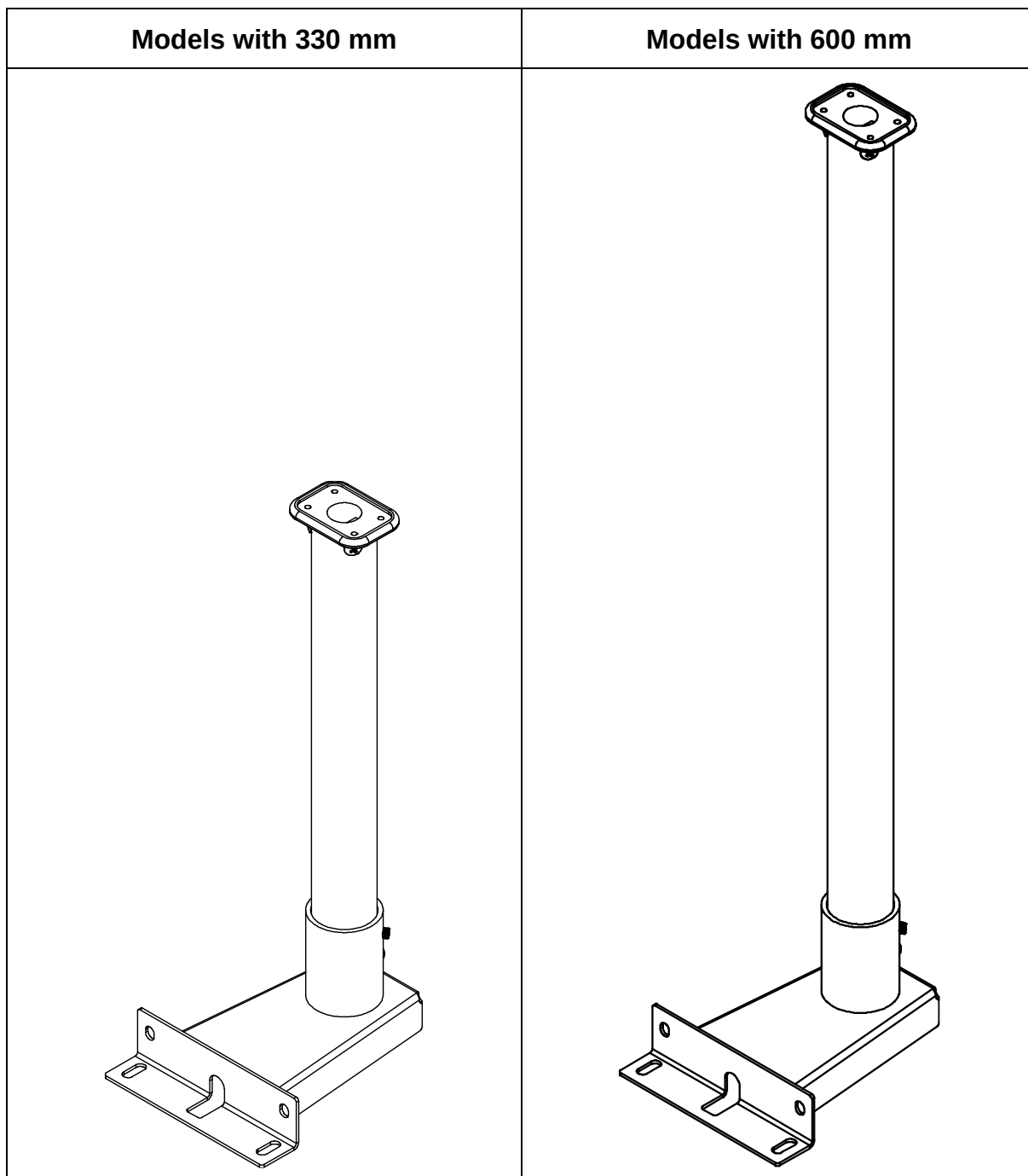







Models with 330 mm:

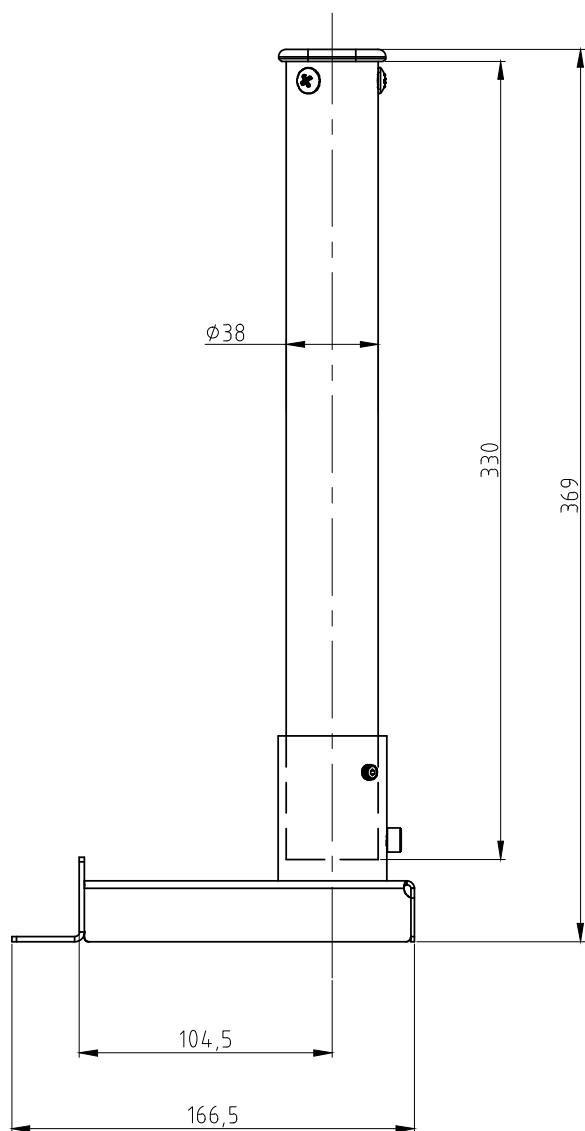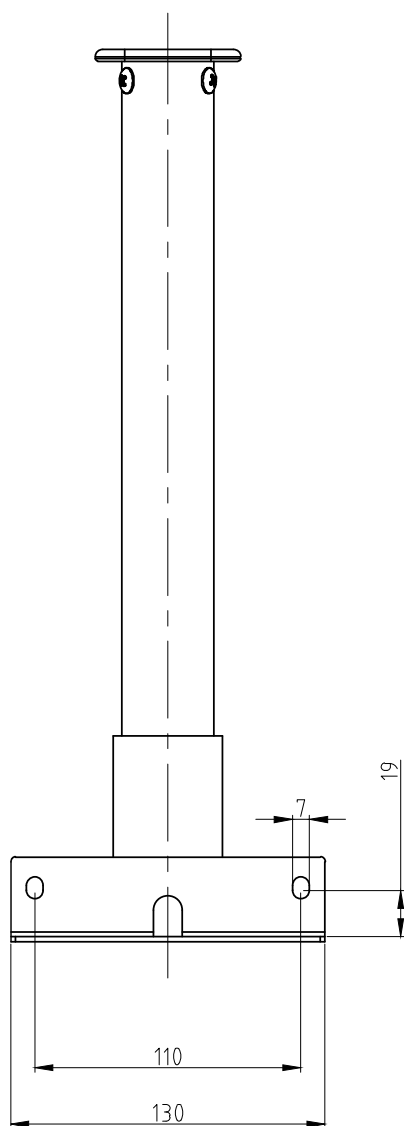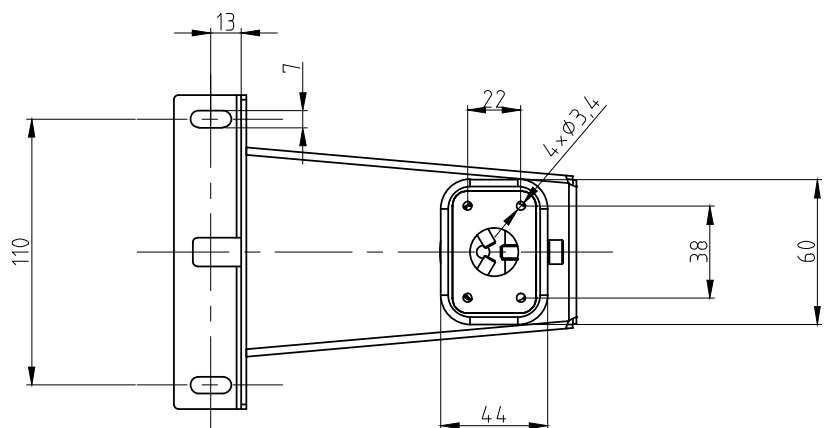

Models with 600 mm:

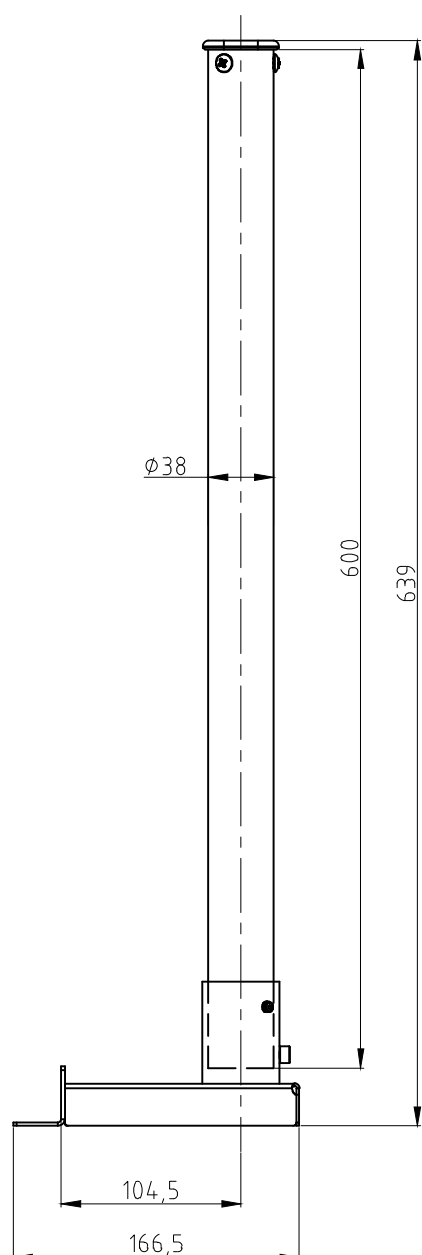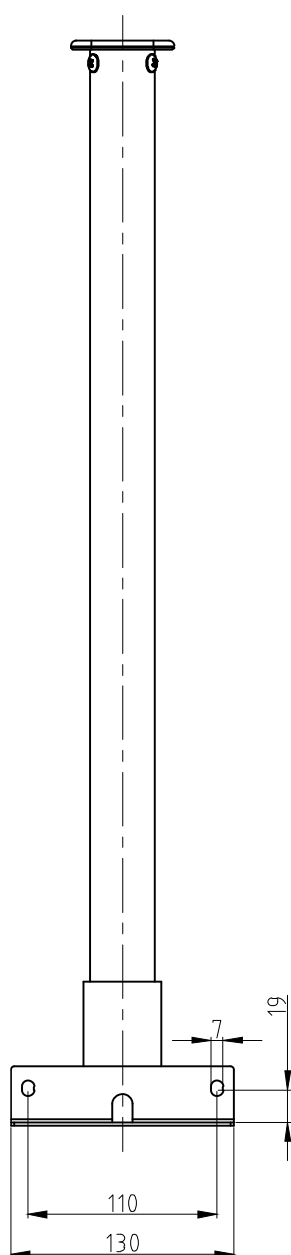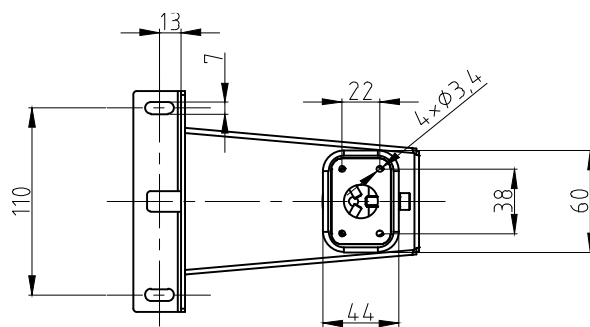

2 Scope of delivery

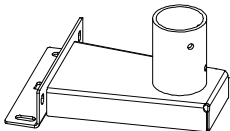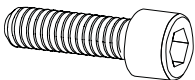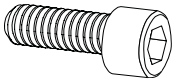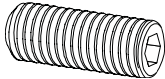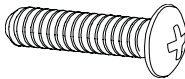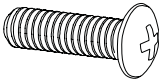
 A Stand foot 1 x	 B Stand tube 1 x	 C Support for display unit 1 x
 S1 M6 x 16 4 x	 S2 Spring washer 4 x	 S3 Washer 4 x
 S4 M6 x 12 1 x	 S5 M6 x 10 2 x	 S6 M5 x 16 3 x
 S7 M4 x 10 4 x		

3 Installation

Assembling the stand foot to the platform:


1. Place the stand foot **(A)** centrally over the perforations on the platform side.
2. First plug the spring washer **(S2)** and then the washer **(S3)** on the screw **(S1)** (4 x).
3. Screw down the stand foot.

Make sure that no cables are jammed or damaged.

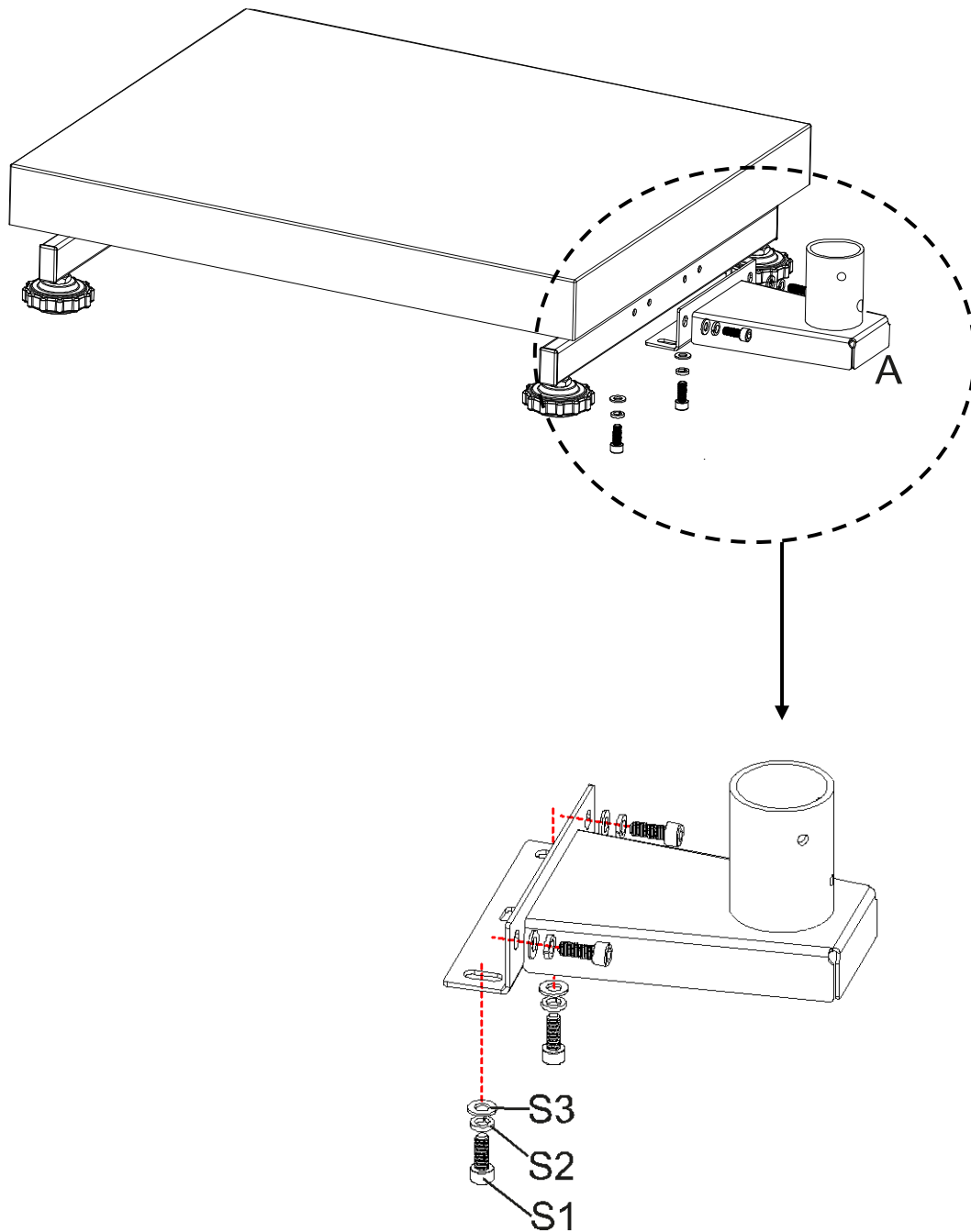

Assembling the stand tube to the stand foot:

4. Unscrew the screws **S4** and **S5** from the stand tube.
5. Push the stand tube **(B)** into the stand foot. Ensure that the perforation of the stand tube matches with the perforation **S4**.
6. Screw down again the screws **S5** and fasten the stand tube using screw **S4**.

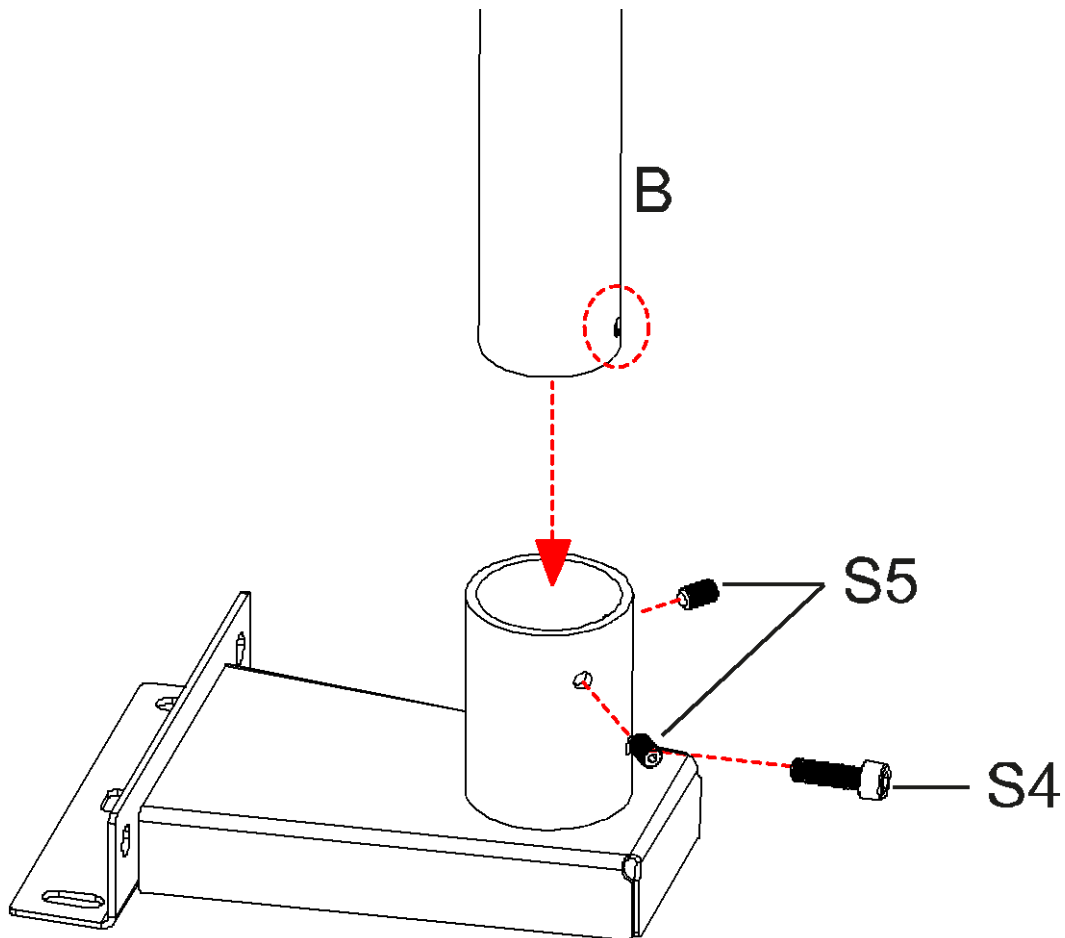

Assembling the display unit to the stand tube:

7. Push the support of the display unit (**C**) into the stand tube and fasten it using screws **S6**.
8. Fasten the display unit to the display unit's support using screws **S7**.

