

PCE INSTRUMENTS

PCE-DHM 5 DIGITAL HANDHELD MICROSCOPE

Deutsch/English/日本語/中文/Español/Français

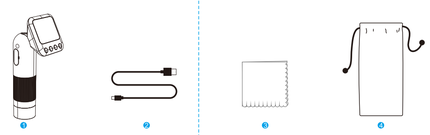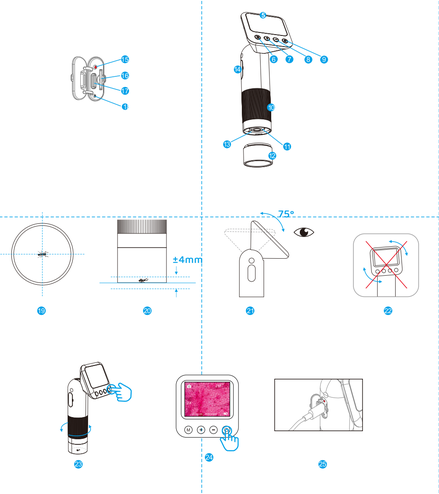Deutsch

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English

Thanks for choosing our microscope. Please read this manual carefully and keep it in a safe place. If you hand the product to others, please be sure to attach this manual. Please heed all warnings and follow all instructions in the manual.

Microscope Kit Introduction

PCE-DHM 5 is a new concept of digital microscope that can be used as a handheld or tabletop magnification tool, to provide bright and clear images, it offers a 6-LED illumination system with adjustable bright white output. To offer flexibility in viewing visual data, the microscope features a tiltable 2" color TFT LCD screen, which can be tilted 75° the most convenient viewing position. and Long-term observation is comfortable for your cervical spine. That can be a great portable companion in your home, outdoor, lab, school, work, or in the field.
Microscope Kit Checklist (This microscope kit includes the following parts)
Microscope and Accessories Use Description

1. The digital Microscope – The microscope holds samples of very small things on clear glass slides. It shines a bright light from its bed plate LED lamp through the small sample,and then lenses make the sample look very big.
2. Type-C data cable - used to transmit data, charge the microscope, and charge the observation platform.
3. Cleaning cloth—only used to wipe some large areas of dirt or water droplets on the microscope body.
4. Carrying case and Straps – This is designed for safe storage of the microscope, shock resistance for transportation, and portability for the viewing experience.

Microscope Details

5. Display - This is where you view and set up and perform microscope shooting functions. The image projected by the acquisition objective lens is displayed on the screen for you to view (you can choose to magnify 1X–5X to view and capture the image/video). At the same time, it also includes the parameter setting of the display and the brightness adjustment of the microscope illumination.
6. M button 7. + button 8. - button 9. Power button

for detailed operation, please refer to the microscope display buttons function and interface description.
10. Focus Ring - Turn this ring very slowly to focus on the specimen you are viewing so you can see it clearly.
11. Objective - It makes the observed object play a role in visual magnification.
12. Transparent plastic cover - mainly for you to clearly see the current state of the specimen when you observe the specimen (the height of the cover is just the best focusing distance of the microscope)
13. LEDs Ring Light - Illuminates the observed object to make it clearly visible.
14. Output Ports - Includes Type-C USB port and TF card slot, allowing you to charge the microscope and store image/video data. (please cover the protective cover when the port is not in use)
15. Microscope charging indicator 16. TF card slot 17. Type-C charging port 18. Reset hole

Microscope Operating Instructions

Pre-observation inspection -Microscope Check
1.1. Long press the⑨on the monitor to light up the screen.
1.2. Please long press the⑦to see if the light is turned on, and short press the⑦/⑧to see if the zoom changes. Press and hold the⑧to see if the light is off.
Begin to observe
2-1. Aim the microscope at the object to be observed, make sure the object is centered on the circular light-transmitting cover⑬and at the edge of the translucent cover⑫.
2-2. Adjust the inclination of the display screen to suit your viewing angle㉑.
The display can be only adjusted at the tilt angle (about 75°) for viewing. please do not rotate it horizontally, otherwise, it may damage the hinge of the display ㉒ .
2-3. Please gently rotate the focus ring of the microscope, and carefully observe the screen changes.If it is clear, you can take pictures or videos㉓㉔.
For more operations, please refer to "Setup Instructions"

Setup Instructions

Microscope display buttons function and interface description
Microscope on/off
Press and hold the power button⑨for 2 seconds, the startup screen will appear. In any mode,
press and hold the power button⑨for 2 seconds to turn off the device.
Photo Mode/Video Modell/Playback Mode
Turn on the device, the default is photo mode. Short press the power button⑨to switch the modes㉘/㉗/㉙.
Photo Mode/Video Modell/Playback Mode
Short press the + button to zoom in, the⑦to zoom out, and the⑨to take pictures/videos. Press and hold the⑦to turn on thelight/brightness+, ⑧brightness-/turn off the light. Long press the⑥to enter the settings menu for this㉘/㉙/㉚.
In playback mode
short press the⑦/⑧to select up and down, and the⑨executes video playback. Press and hold the⑦/⑧to fast forward/rewind during playback. Long press the⑨to delete the file.
After entering the setting menu, short press the⑥to switch the setting type (photo/video setting or general setting), the⑦/⑧can be selected up and down, and the power button is executed to confirm. When you have enteredthe menu bar, short press the⑥to return to the previous menu.
㉖ Photography interface ㉗ Video interface ㉘ Backflow interface
㉙ Photo settings interface ㉚ Video settings interface ㉛ General settings interface

Other instructions

Charging Instructions
Microscope charging: Insert the type-c plug into the Type-c port of the microscope, and connect the other end to the USB power supply for charging ㉕ The charging time is about 2 hours, and the usage time is about 2 hours. (Red light when charging, green light when fully charged).
Frequently Asked Questions

Problem	Cause	Solution
Microscope does not turn on	The battery may be dead	Please charge.
The picture is noisy and dark	No lights on	Please turn on the light
system stuck	Press the reset key	Please insert a needle into the reset hole and press it to reset the system.
The picture is not clear	The focus is not in place, the sample is notplaced or handled properly, and the surface of the slide is not clean.	Please check whether the focus is correct, whether the specimen is placed properly, and whether the slide is clean and free of foreign matter.
The screen is too bright	The exposure parameters are not adjusted correctly.	Please enter the setting menu to adjust exposure parameters

Precautions
• Keep children and animals away from any experiments or programs.
• Store the microscope unit out of the reach of children.
• WARNING: This microscope contains glass parts and is intended for use by adults only.
• This product is a precision instrument, please avoid dropping, bumping or receiving heavy blows.
• Do not place this product in a place where the temperature is too high, the high temperature may damage or explode the battery inside the product.

Friendly Reminder

• Please clean this product with a soft, dry cloth.
• Contact authorized service personnel when service is required.
• This product contains optical components and precision electronic circuits, which may be disassembled without permission, resulting in failure.
• Damage to the machine due to improper use is not covered by the warranty.

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